



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

# PRODUCT INFORMATION



## NAADP (sodium salt)

Item No. 25992

**Formal Name:** P'→5'-ester with 3-carboxy-1-β-D-ribofuranosylpyridinium, 2'-(dihydrogen phosphate), adenosine 5'-(trihydrogen diphosphate), inner salt, tetrasodium salt  
**Synonym:** Nicotinic acid adenine dinucleotide phosphate


**MF:** C<sub>21</sub>H<sub>23</sub>N<sub>6</sub>O<sub>18</sub>P<sub>3</sub> • 4Na  
**FW:** 832.3

**Purity:** ≥95%

**Supplied as:** A solid

**Storage:** -20°C

**Stability:** ≥2 years

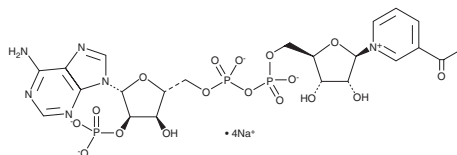

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

### Laboratory Procedures

NAADP (sodium salt) is supplied as a solid. A stock solution may be made by dissolving the NAADP (sodium salt) in water. The solubility of NAADP (sodium salt) in water is approximately 50 mg/ml. We do not recommend storing the aqueous solution for more than one day.

### Description

Nicotinic acid adenine dinucleotide phosphate (NAADP) is a secondary messenger that induces calcium mobilization.<sup>1,2</sup> It induces calcium release from endosomes and lysosomes *via* two-pore channel 2 (TPC2) and TPC1, which then stimulates large-scale calcium release from granules and the endoplasmic reticulum mediated by type 1 ryanodine receptors (RyR2s), RyR3s, and inositol-(1,4,5)-triphosphate receptors (IP<sub>3</sub>Rs).<sup>1,3</sup> NAADP induces calcium mobilization in sea urchin and starfish eggs post fertilization to block polyspermy and activate embryogenesis.<sup>3</sup> NAADP-induced calcium mobilization induces VEGF-mediated angiogenesis in human umbilical vein endothelial cells (HUVECs). It also alkalizes lysosomal pH thereby inhibiting fusion between autophagosomes and lysosomes and arresting autophagic flux in mouse embryonic stem cells.<sup>4</sup>

### References

1. Gerasimenko, J.V., Charlesworth, R.M., Sherwood, M.W., *et al.* Both RyRs and TPCs are required for NAADP-induced intracellular Ca<sup>2+</sup> release. *Cell Calcium* **58(3)**, 237-245 (2015).
2. Lee, H.C. Mechanisms of calcium signaling by cyclic ADP-ribose and NAADP. *Physiological Reviews* **77(4)**, 1133-1164 (1997).
3. Parrington, J., Lear, P., and Hachem, A. Calcium signals regulated by NAADP and two-pore channels--their role in development, differentiation and cancer. *Int. J. Dev. Biol.* **59(7-9)**, 341-355 (2015).
4. Lu, Y., Hao, B., Graeff, R., *et al.* NAADP/TPC2/Ca<sup>2+</sup> signaling inhibits autophagy. *Commun. Integr. Biol.* **6(6)**, e27595 (2013).

#### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

#### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

#### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 09/21/2018

#### CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA

**PHONE:** [800] 364-9897  
[734] 971-3335

**FAX:** [734] 971-3640

CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM