



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



Dimethyldioctadecylammonium (bromide)

Item No. 26120

CAS Registry No.: 3700-67-2
Formal Name: N,N-dimethyl-N-octadecyl-1-octadecanaminium, monobromide
Synonyms: DDA, DODAB
MF: C₃₈H₈₀N • Br
FW: 631.0
Purity: ≥95%
Supplied as: A crystalline solid
Storage: -20°C
Stability: ≥2 years

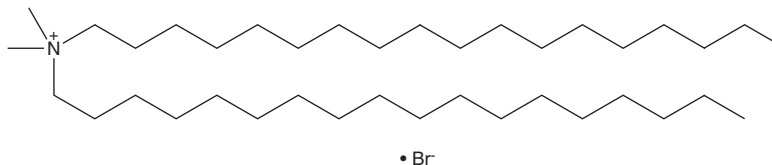

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Dimethyldioctadecylammonium (DDA) (bromide) is supplied as a crystalline solid. A stock solution may be made by dissolving the DDA (bromide) in the solvent of choice. DDA (bromide) is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide, which should be purged with an inert gas. The solubility of DDA (bromide) in these solvents is approximately 33, 16, and 2 mg/ml, respectively.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. Organic solvent-free aqueous solutions of DDA (bromide) can be prepared by directly dissolving the crystalline solid in aqueous buffers. The solubility of DDA (bromide) in PBS, pH 7.2, is approximately 2 mg/ml. We do not recommend storing the aqueous solution for more than one day.

Description

DDA is a cationic amphipathic lipid.^{1,2} DDA liposomes containing an Ag85B-ESAT-6 antigen induce antigen deposition at an intramuscular or subcutaneous injection site in mice, increasing immune cell exposure to the antigen.³ In a guinea pig model of *M. tuberculosis* infection, spleen bacterial load is reduced and lung and spleen lesion numbers are decreased when the mycobacterial lipid antigens Ac₂SGL and PIM₂ are administered in liposomes comprised of DDA and a synthetic analog of the mycobacterial cord factor trehalose 6,6-dibehenate (TDB).² DDA has also been used in the study of lipid bilayer dynamics.⁴

References

1. De Serrano, L.O. and Burkhart, D.J. Liposomal vaccine formulations as prophylactic agents: Design considerations for modern vaccines. *J. Nanobiotechnology* **15**(1), 83 (2017).
2. Larrouy-Maumus, G., Layre, F., Clark, S., et al. Protective efficacy of a lipid antigen vaccine in a guinea pig model of tuberculosis. *Vaccine* **35**(10), 1395-1402 (2017).
3. Henriksen-Lacey, M., Bramwell, V.W., Christensen, D., et al. Liposomes based on dimethyldioctadecylammonium promote a depot effect and enhance immunogenicity of soluble antigen. *J. Control Release* **142**(2), 180-186 (2010).
4. Dubey, P.S., Sharma, V.K., Srinivasan, H., et al. Effects of NSAIDs on the dynamics and phase behavior of DODAB bilayers. *J. Phys. Chem. B.* **122**(43), 9962-9972 (2018).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/13/2018

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD

ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM