



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com


[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)

PRODUCT INFORMATION



Dofetilide-d₄ Item No. 26451

CAS Registry No.: 1189700-56-8
Formal Name: N-[4-[2-[methyl[2-[4-[(methylsulfonyl)amino]phenoxy]ethyl-1,1,2,2-d₄]amino]ethyl]phenyl]-methanesulfonamide
MF: C₁₉H₂₃D₄N₃O₅S₂
FW: 445.6
Chemical Purity: ≥98% (Dofetilide)
Deuterium Incorporation: ≥99% deuterated forms (d₁-d₄); ≤1% d₀
Supplied as: A solid
Storage: -20°C
Stability: ≥2 years

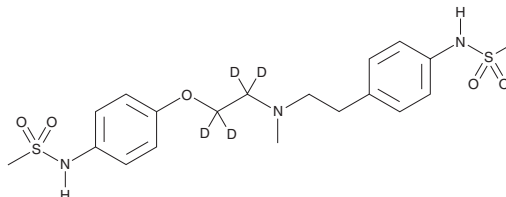

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Dofetilide-d₄ is intended for use as an internal standard for the quantification of dofetilide (Item No. 15045) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

Dofetilide-d₄ is supplied as a solid. A stock solution may be made by dissolving the dofetilide-d₄ in the solvent of choice, which should be purged with an inert gas. Dofetilide-d₄ is slightly soluble in chloroform and methanol.

Description

Dofetilide is a class III antiarrhythmic agent that prolongs cardiac action potential duration by selectively inhibiting the rapidly activating inward rectifying component of net delayed rectifier K⁺ current (IC₅₀ = 31.5 nM in guinea pig cardiomyocytes).¹ However, at 1 μM, dofetilide has pro-arrhythmic activity, inducing early afterdepolarizations (prolonged repolarization) in cell models and Torsade de Pointes in a rabbit screen for proarrhythmic properties when administered at a dose of 10 mg/kg.² Formulations containing dofetilide have been used in the treatment of highly symptomatic atrial fibrillation and the conversion of flutter to normal sinus rhythm.³

References

1. Jurkiewicz, N.K. and Sanguinetti, M.C. Rate-dependent prolongation of cardiac action potentials by a methanesulfonamide class III antiarrhythmic agent. *Circ. Res.* **72(1)**, 75-83 (1993).
2. Nalos, L., Varkevisser, R., Jonsson, M.K., *et al.* Comparison of the I_{Kr} blockers moxifloxacin, dofetilide and E-4031 in five screening models of pro-arrhythmia reveals lack of specificity of isolated cardiomyocytes. *Br. J. Pharmacol.* **165(2)**, 467-478 (2012).
3. Dunnink, A., van Opstal, J.M., Oosterhoff, P., *et al.* Ventricular remodelling is a prerequisite for the induction of dofetilide-induced torsade de pointes arrhythmias in the anaesthetized, complete atrio-ventricular-block dog. *Europace* **14(3)**, 431-436 (2012).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/30/2018

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD

ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM