

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



Troglitazone-d₄

Item No. 28908

Formal Name: 5-[[4-[(3,4-dihydro-6-hydroxy-2,5,7,8-

tetramethyl-2H-1-benzopyran-2-yl)methoxyl

phenyl-d₄]methyl]-2,4-thiazolidinedione

MF: $C_{24}H_{23}D_4NO_5S$

FW: 445.6


Chemical Purity: ≥95% (Troglitazone)

Deuterium

Incorporation: ≥99% deuterated forms (d_1-d_4) ; ≤1% d_0

Supplied as: A solid Storage: -20°C Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

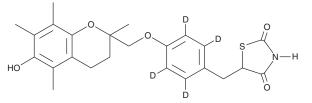

Laboratory Procedures

Troglitazone- d_{Δ} is intended for use as an internal standard for the quantification of troglitazone (Item No. 71750) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

Troglitazone- d_A is supplied as a solid. A stock solution may be made by dissolving the troglitazone- d_A in the solvent of choice, which should be purged with an inert gas. Troglitazone-d₁ is soluble in methanol.

Description

Troglitazone is a selective agonist of peroxisome proliferator-activated receptor γ (PPARγ; EC₅₀s = 0.78 and 0.55 μ M for the mouse and human receptors, respectively, in a transactivation assay). ¹ It is selective for PPAR γ over PPAR α and PPAR δ , at which it is inactive at concentrations up to 10 μ M. Troglitazone (500 mg/kg twice per day) exhibits antihyperglycemic and antihyperlipidemic properties in Zucker diabetic fatty rats, reducing plasma glucose and triglyceride levels by 61 and 87%, respectively.² It also induces cell cycle arrest at the G₁ phase in SK-HEP-1 and Hep3B hepatocellular carcinoma cells when used at a concentration of 10 µM and increases apoptosis in these cells at concentrations of 30 µM and higher.3

References

- 1. Willson, T.M., Brown, P.J., Sternbach, D.D., et al. The PPARs: From orphan receptors to drug discovery. J. Med. Chem. 43(4), 527-550 (2000).
- 2. Henke, B.R., Blanchard, S.G., Brackeen, M.F., et al. N-(2-Benzoylphenyl)-L-tyrosine PPARy agonists. 1. Discovery of a novel series of potent antihyperglycemic and antihyperlipidemic agents. J. Med. Chem. 41(25), 5020-5036 (1998).
- 3. Yoshizawa, K., Cioca, D.P., Kawa, S., et al. Peroxisome proliferator-activated receptor y ligand troglitazone induces cell cycle arrest and apoptosis of hepatocellular carcinoma cell lines. Cancer 95(10), 2243-2251 (2002).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information Buyer agrees to purchase the material can be found on our website.

Copyright Cayman Chemical Company, 10/22/2019

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM