



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



(±)5(6)-EET-d₁₁ methyl ester Item No. 29382

Formal Name: methyl 4-((2S,3R)-3-((2Z,5Z,8Z)-tetradeca-2,5,8-trien-1-yl-oxiran-2-yl)butanoate

Synonym: (±)5,6-EpETrE-d₁₁ methyl ester

MF: C₂₁H₂₃D₁₁O₃

FW: 345.6


Chemical Purity: ≥98% ((±)5(6)-EET methyl ester)

Deuterium Incorporation: ≥99% deuterated forms (d₁-d₁₁); ≤1% d₀

Supplied as: A solution in ethanol

Storage: -20°C

Stability: ≥1 year

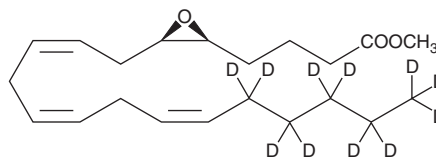

NOTE: Relative stereochemistry shown in chemical structure

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

(±)5(6)-EET-d₁₁ methyl ester is intended for use as an internal standard for the quantification of 5(6)-EET methyl ester by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

(±)5(6)-EET-d₁₁ methyl ester is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. The solubility of (±)5(6)-EET-d₁₁ methyl ester in these solvents is approximately 50 mg/ml.

Description

(±)5(6)-EET methyl ester (Item No. 50210) is an esterified form of (±)5(6)-EET (Item No. 50211). 5(S),6(R)-EET and 5(R),6(S)-EET are formed *via* epoxidation of arachidonic acid (Item Nos. 90010 | 90010.1 | 10006607) by a variety of cytochrome P450 (CYP) isoforms, including CYP2B1, CYP2B2, CYP2C11, CYP2C23, and CYP2C24.^{1,2}

References

1. Spector, A.A. and Norris, A.W. Action of epoxyeicosatrienoic acids on cellular function. *Am. J. Physiol. Cell Physiol.* **292**(3), C996-C1012 (2007).
2. Capdevila, J.H., Falck, J.R., and Harris, R.C. Cytochrome P450 and arachidonic acid bioactivation: Molecular and functional properties of the arachidonate monooxygenase. *J. Lipid Res.* **41**(2), 163-181 (2000).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 02/26/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM