



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



5-(N,N-hexamethylene)-Amiloride

Item No. 29788

CAS Registry No.: 1428-95-1
Formal Name: 3-amino-N-(aminoiminomethyl)-6-chloro-5-(hexahydro-1H-azepin-1-yl)-2-pyrazinecarboxamide
Synonym: HMA
MF: C₁₂H₁₈ClN₇O
FW: 311.8
Purity: ≥98%
UV/Vis.: λ_{max}: 232, 295, 375 nm
Supplied as: A crystalline solid
Storage: -20°C
Stability: ≥2 years

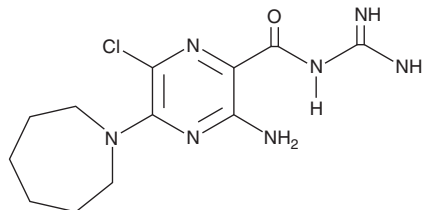

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

5-(N,N-hexamethylene)-Amiloride (HMA) is supplied as a crystalline solid. A stock solution may be made by dissolving the HMA in the solvent of choice, which should be purged with an inert gas. HMA is soluble in organic solvents such as DMSO and dimethyl formamide. The solubility of HMA in these solvents is approximately 10 and 3 mg/ml, respectively.

HMA is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, HMA should first be dissolved in DMSO and then diluted with the aqueous buffer of choice. HMA has a solubility of approximately 0.2 mg/ml in a 1:4 solution of DMSO:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

HMA is a derivative of amiloride (Item No. 14409) with diverse biological activities.¹⁻⁵ It is an allosteric antagonist of adenosine A_{2A} receptors (K_i = 3.3 μM).² HMA inhibits the cation-selective ion channel formed by the HIV-1 viral protein Vpu when used at a concentration of 50 μM, as well as budding of virus-like particles in HeLa cells expressing the HIV-1 proteins Gag and Vpu when used at a concentration of 10 μM.¹ It also blocks the cation-selective ion channels formed by the hepatitis C virus (HCV) protein p7.³ HMA (40 μM) induces necrosis in and reduces the viability of MCF-7, MDA-MB-231, T47D, SK-BR-3, Met-1, and NDL breast cancer cells but not cardiomyocytes or uterine, pulmonary, and renal epithelial cells.⁴ HMA protects against post-ischemic contractile dysfunction and reduces coronary effluent creatine phosphokinase activity in a model of ischemia-reperfusion injury using isolated rat right ventricular free walls.⁵

References

1. Ewart, G.D., Mills, K., Cox, G.B., et al. *Eur. Biophys. J.* **31(1)**, 26-35 (2002).
2. Gao, Z.-G. and Ijzerman, A.P. *Biochem. Pharmacol.* **60(5)**, 669-679 (2000).
3. Premkumar, A., Wilson, L., Ewart, G.D., et al. *FEBS Lett.* **557(1-3)**, 99-103 (2004).
4. Rowson-Hodel, A.R., Berg, A.L., Wald, J.H., et al. *Cancer Lett.* **375(1)**, 62-72 (2016).
5. Meng, H.-P., Maddaford, T.G., and Pierce, G.N. *Am. J. Physiol.* **264(6 Pt. 2)**, H1831-H1835 (1993).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 02/26/2020

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM