



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



Oxaloacetic Acid

Item No. 30280

CAS Registry No.: 328-42-7
Formal Name: 2-oxo-butanedioic acid
Synonyms: NSC 284205, NSC 77688, OAA, Oxalacetic Acid
MF: C₄H₄O₅
FW: 132.1
Purity: ≥95%
UV/Vis.: λ_{max}: 259 nm
Supplied as: A solid
Storage: -20°C
Stability: ≥4 years

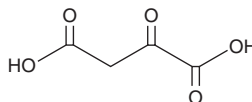

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Oxaloacetic acid is supplied as a solid. A stock solution may be made by dissolving the oxaloacetic acid in the solvent of choice, which should be purged with an inert gas. Oxaloacetic acid is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide. The solubility of oxaloacetic acid in these solvents is approximately 15 mg/ml.

Description

Oxaloacetic acid is an α-keto acid and a key component of cellular metabolism in its conjugate base form, oxaloacetate.¹⁻³ Oxaloacetate reacts with acetyl-coenzyme A (acetyl-CoA; Item No. 16160) and water to form citrate in the first step of the citric acid cycle and is regenerated by oxidation of L-malate in the final step.¹ It is an intermediate in gluconeogenesis that is formed in mitochondria via carboxylation of pyruvate and subsequently decarboxylated and phosphorylated to form phosphoenolpyruvate.² It can be converted to aspartate via addition of an amino group from glutamate.³ Oxaloacetate (30 μmol/min per 100 g for 30 minutes, i.v.) reduces blood glutamate levels, severity of neurological dysfunction, and brain edema in a rat model of closed head injury.⁴

References

1. Berg, J.M., Tymoczko, J.L., and Stryer, L. The citric acid cycle oxidizes two-carbon units. *Biochemistry*, 5th edition, W H Freeman (2002).
2. Berg, J.M., Tymoczko, J.L., and Stryer, L. Glucose can be synthesized from noncarbohydrate precursors. *Biochemistry*, 5th edition, W H Freeman (2002).
3. Berg, J.M., Tymoczko, J.L., and Stryer, L. Amino acids are made from intermediates of the citric acid cycle and other major pathways. *Biochemistry*, 5th edition, W H Freeman (2002).
4. Zlotnik, A., Gurevich, B., Tkachov, S., et al. Brain neuroprotection by scavenging blood glutamate. *Exp. Neurol.* 203(1), 213-220 (2007).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 12/08/2022

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD

ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM