

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



Norgestimate-d₄

Item No. 31176

CAS Registry No.: 1263194-12-2

Formal Name: $(17\alpha)-17$ -(acetyloxy)-13-ethyl-18,19-

dinorpregn-4-en-20-yn-3-one-

2,2,4,6,6,10-d₆, 3-oxime

Synonym: Dexnorgestrel acetime-d₆

C₂₃H₂₅D₆NO₃ 375.5 MF:

FW:


Chemical Purity: ≥95% (mixture of isomers) (Norgestimate)

Deuterium

Incorporation: \geq 99% deuterated forms (d₁-d₆); \leq 1% d₀

Supplied as: A solid -20°C Storage: Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

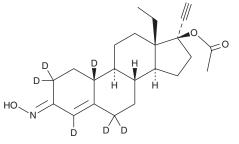

Laboratory Procedures

Norgestimate-d₆ is intended for use as an internal standard for the quantification of norgestimate (Item No. 16547) by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

Norgestimate- d_6 is supplied as a solid. A stock solution may be made by dissolving the norgestimate- d_6 in the solvent of choice, which should be purged with an inert gas. Norgestimate-d, is soluble in the organic solvent methanol.

Description

Norgestimate is a progestin and a prodrug form of norelgestromin and levonorgestrel (Item No. 10006318).^{1,2} It binds to progestin receptors (mean IC_{50} = 3.48 nM for rabbit uterine receptors) and is selective for progestin over androgen receptors (mean IC_{50} = 764 nM for rat prostatic receptors) in radioligand binding assays.³ Norgestimate maintains pregnancy in ovariectomized female rats in a dose-dependent manner. 1 It inhibits ovulation in 100% of female rats when administered at a dose of 0.5 mg/kg. Norgestimate (2 mg/kg) completely inhibits estrone-induced vaginal cornification in ovariectomized female rats, indicating antiestrogenic activity. Formulations containing norgestimate in combination with ethinyl estradiol have been used as oral contraceptives and in the treatment of acne vulgaris.

References

- 1. Hahn, D.W., Allen, G.O., and McGuire, J.L. The pharmacological profile of norgestimate, a new orally active progestin. Contraception 16(5), 541-553 (1977).
- 2. Stanczyk, F.Z. All progestins are not created equal. Steroids 68(10-13), 879-890 (2003).
- 3. Phillips, A., Demarest, K., Hahn, D.W., et al. Progestational and androgenic receptor binding affinities and in vivo activities of norgestimate and other progestins. Contraception 41(4), 399-410 (1990).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information Buyer agrees to purchase the material can be found on our website.

Copyright Cayman Chemical Company, 08/06/2020

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM