

# Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



# Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

# Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

## SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



# PRODUCT INFORMATION



12(S)-HETE-d<sub>s</sub> Item No. 334570

CAS Registry No.: 84807-90-9

Formal Name: 12S-hydroxy-5Z,8Z,10E,14Z-

eicosatetraenoic-5,6,8,9,11,12,14,15- $d_8$  acid

MF:

FW: 328.5 **Chemical Purity:** ≥98% (12(S)-HETE)

Deuterium


Incorporation: ≥99% deuterated forms  $(d_1-d_8)$ ; ≤1%  $d_0$ 

UV/Vis.:  $\lambda_{max}$ : 237 nm

Supplied as: A solution in acetonitrile

-20°C Storage: Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

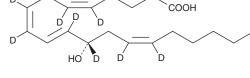

#### **Laboratory Procedures**

12(S)-HETE-d<sub>8</sub> is intended for use as an internal standard for the quantification of 12-HETE by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

12(S)-HETE-d<sub>8</sub> is supplied as a solution in acetonitrile. To change the solvent, simply evaporate the acetonitrile under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as ethanol, DMSO, and dimethyl formamide purged with an inert gas can be used. 12(S)-HETE-do is miscible in these solvents.

## Description

(±)12-HETE is formed via non-enzymatic oxidation of arachidonic acid (Item Nos. 90010 | 90010.1 10006607).1,2 12(S)- and 12(R)-HETE are formed by 12(S)- and 12(R)-lipoxygenase-mediated oxidation of arachidonic acid, respectively.<sup>3,4</sup> 12(R)-HETE can also be formed by oxidation of arachidonic acid mediated by cytochrome P450s (CYP450s).1

#### References

- 1. Powell, W.S. and Rokach, J. Biosynthesis, biological effects, and receptors of hydroxyeicosatetraenoic acids (HETEs) and oxoeicosatetraenoic acids (oxo-ETEs) derived from arachidonic acid. Biochim. Biophys. Acta 1851(4), 340-355 (2014).
- 2. Derogis, P.B.M.C., Chaves-Filho, A.B., and Miyamoto, S. Characterization of hydroxy and hydroperoxy polyunsaturated fatty acids by mass spectrometry. Bioactive lipids in health and disease. Trostchansky, A. and Rubbo, H., editors, Springer (2019).
- Kühn, H., Banthiya, S., and van Leyen, K. Mammalian lipoxygenases and their biological relevance. Biochim. Biophys. Acta 1851(4), 308-330 (2015).
- 4. Bürger, F., Krieg, P., Marks, F., et al. Positional- and stereo-selectivity of fatty acid oxygenation catalysed by mouse (12S)-lipoxygenase isoenzymes. Biochem J. 348(Pt 2), 329-335 (2000).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

#### WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information Buyer agrees to purchase the material can be found on our website.

Copyright Cayman Chemical Company, 03/17/2021

### **CAYMAN CHEMICAL**

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM