



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



9(10)-Nitrooleate

Item No. 33896

CAS Registry No.: 1092676-99-7
Formal Name: 9-nitro-9E-octadecenoic acid, 10-nitro-9E-octadecenoate
Synonyms: 9(10)-Nitrooleic Acid, OA-NO₂, 9(10)-nitro-9-*trans*-Octadecenoic Acid
MF: C₁₈H₃₃NO₄ • C₁₈H₃₃NO₄
FW: 654.9
Purity: ≥95% (1:1 mixture of isomers)
Supplied as: A solution in ethanol
Storage: -20°C
Stability: ≥1 year

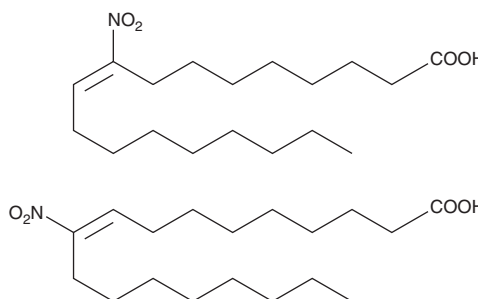

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Description

9(10)-Nitrooleate is a mixture of the nitroalkene endogenous lipid signaling molecules 9-nitrooleate (Item No. 10008042) and 10-nitrooleate (Item No. 10008043). 9(10)-Nitrooleate is formed by nitration of oleic acid (Item Nos. 90260 | 24659) by peroxyxynitrite, acidified nitrite, and myeloperoxidase in the presence of hydrogen peroxide (H₂O₂) and nitrite.¹ It induces transactivation of peroxisome proliferator-activated receptor γ (PPAR γ) or nuclear factor erythroid 2-related factor 2 (Nrf2) in reporter assays using MCF-7 cells when used at concentrations greater than or equal to 0.3 and greater than or equal to 3 μ M, respectively.² 9(10)-Nitrooleate also induces transactivation of PPAR γ , PPAR δ , or PPAR α in reporter assays using CV-1 cells in a concentration-dependent manner.¹ It induces adipogenesis in 3T3-L1 preadipocytes and glucose uptake in differentiated 3T3-L1 adipocytes when used at a concentration of 3 μ M. *In vivo*, 9(10)-nitrooleate (500 μ g/kg every six hours for 24 hours) decreases ischemia-induced increases in plasma levels of blood urea nitrogen, creatinine, and TNF- α and kidney levels of myeloperoxidase (Mpo) and IL-1 β and reduces tubular necrosis and tubule dilation in a mouse model of renal ischemia-reperfusion injury.³

References

1. Baker, P.R., Lin, Y., Schopfer, F.J., *et al.* Fatty acid transduction of nitric oxide signaling. Multiple nitrated unsaturated acid derivatives exist in human blood and urine serve as endogenous peroxisome proliferator-activated receptor ligands. *J. Biol. Chem.* **280**(51), 42464-42475 (2005).
2. Bates, D.J.P., Smitherman, P.K., Townsend, A.L., *et al.* Nitroalkene fatty acids mediate activation of Nrf2/ARE-dependent and PPAR γ -dependent transcription by distinct signaling pathways and with significantly different potencies. *Biochemistry* **50**(36), 7765-7773 (2011).
3. Liu, H., Jia, Z., Soodvilai, S., *et al.* Nitro-oleic acid protects the mouse kidney from ischemia and reperfusion injury. *Am. J. Physiol. Renal. Physiol.* **295**(4), F942-F949 (2008).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 10/01/2024

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM