

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


Product Information



5(R)-HETE

Item No. 34225

CAS Registry No.:	61641-47-2
Formal Name:	5R-hydroxy-6E,8Z,11Z,14Z-
	eicosatetraenoic acid
MF:	C ₂₀ H ₃₂ O ₃
FW:	320.5
Purity:	≥98%
Stability:	≥1 year at -20°C
Supplied as:	A solution in ethanol
UV/Vis:	λ _{max} : 236 nm ε: 27,000
Misc:	Oxygen and light sensitive

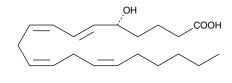

Laboratory Procedures

For long term storage, we suggest that 5(R)-HETE be stored as supplied at -20°C. It should be stable for at least one year.

5(R)-HETE is supplied as a solution in ethanol. To change the solvent, evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO or dimethyl formamide purged with an inert gas can be used. The solubility of 5(R)-HETE in these solvents is approximately 50 mg/ml.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made just prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. If an organic solvent-free solution of 5(R)-HETE is needed, evaporate the ethanol under a gentle stream of nitrogen and dissolve the neat oil in the buffer of choice. The solubility of 5(R)-HETE in PBS, pH 7.2, is approximately 1 mg/ml. Store aqueous solutions of 5(R)-HETE on ice and use within twelve hours.

5(R)-HETE is a monohydroxy polyunsaturated fatty acid produced by the nonenzymatic oxidation of arachidonic acid. It is a potent chemotactic agent. 5(R)-HETE induces human neutrophil migration across cellular and noncellular barriers in a dose- and time-dependent fashion, whereas 5(S)-HETE is inactive.¹

Reference

1. Bittleman, D.B. and Casale, T.B. 5-Hydroxyeicosatetraenoic acid (HETE)-induced neutrophil transcellular migration is dependent upon enantiomeric structure. Am. J. Respir. Cell Mol. Biol. 12, 260-267 (1995).

Cayman Chemical

Mailing address

1180 E. Ellsworth Road Ann Arbor, MI 48108 USA

Phone (800) 364-9897 (734) 971-3335

Fax (734) 971-3640

E-Mail

custserv@caymanchem.com

Web www.caymanchem.com

WARNING: This product is for laboratory research only: not for administration to humans. Not for human or veterinary DIAGNOSTIC OR THERAPEUTIC USE.

MATERIAL SAFETY DATA

This material should be considered hazardous until information to the contrary becomes available. Do not ingest, swallow, or inhale. Do not get in eyes, on skin, or on clothing. Wash thoroughly after handling. This information contains some, but not all, of the information required for the safe and proper use of this material. Before use, the user must review the complete Material Safety Data Sheet, which has been sent via email to your institution

WARRANTY AND LIMITATION OF REMEDY

Cayman Chemical Company makes no warranty or guarantee of any kind, whether written or oral, expressed or implied, including without limitation, any warranty of fitness for a particular purpose, suitability and merchantability, which extends beyond the description of the chemicals hereof. Cayman warrants only to the original customer that the material will meet our specifications

purpose, suitability and merchanicality, which exclude ocycles are construction of the time of delivery. at the time of delivery. Cayman will carry out its delivery obligations with due care and skill. Thus, in no event will Cayman have **any obligation or liability**, whether in tort (including negligence) or in contract, for any direct, indirect, incidental or consequential damages, even if Cayman is informed about their possible existence. This limitation of liability does not apply in the case of intentional acts or negligence of Cayman, its directors or its employees. Buyer's **exclusive remedy** and Cayman's sole liability hereunder shall be limited to a <u>refund</u> of the purchase price, or at Cayman's option, the <u>replacement</u>, at no cost to Buyer, of all material that the purchase price of the purchase price of the purchase price.

Said refund or replacement is conditioned on Buyer giving written notice to Cayman within thirty (30) days after arrival of the material at its destination. Failure of Buyer to give said notice within try (30) days shall constitute a waiver by Buyer of all claims hereunder with respect to said material. For further details, please refer to our Warranty and Limitation of Remedy located on our website and in our catalog.