



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


# PRODUCT INFORMATION



## 5-OxoETE

Item No. 34250

**CAS Registry No.:** 106154-18-1  
**Formal Name:** 5-oxo-6E,8Z,11Z,14Z-eicosatetraenoic acid  
**Synonym:** 5-KETE  
**MF:** C<sub>20</sub>H<sub>30</sub>O<sub>3</sub>  
**FW:** 318.5  
**Purity:** ≥95%  
**UV/Vis.:** λ<sub>max</sub>: 279 nm  
**Supplied as:** A solution in ethanol  
**Storage:** -80°C  
**Stability:** ≥1 year

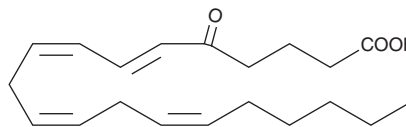

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

### Laboratory Procedures

5-oxoETE is supplied as a solution in ethanol. To change the solvent, evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. 5-oxoETE is miscible in these solvents.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. If an organic solvent-free solution of 5-oxoETE is needed, it can be prepared by evaporating the ethanol and directly dissolving the neat oil in aqueous buffers. The solubility of 5-oxoETE in PBS, pH 7.2, is approximately 0.8 mg/ml. For greater aqueous solubility, 5-oxoETE can be directly dissolved in 0.1 M Na<sub>2</sub>CO<sub>3</sub> (solubility of 2 mg/ml) and then diluted with PBS (pH 7.2) to achieve the desired concentration or pH. Store aqueous solutions of 5-oxoETE on ice and use within 12 hours of preparation. We do not recommend storing the aqueous solution for more than one day.

### Description

5-OxoETE is a polyunsaturated keto acid formed by the oxidation of 5-HETE in human neutrophils by a specific dehydrogenase.<sup>1</sup> It stimulates cytosolic calcium levels in neutrophils with an EC<sub>50</sub> value of 2 nM.<sup>2</sup> 5-OxoETE selectively stimulates the migration and degranulation of eosinophils and activates the MAPK pathway in stimulated neutrophils *via* a specific G protein-coupled receptor.<sup>3-6</sup>

### References

1. Powell, W.S., Gravelle, F., and Gravel, S. *J. Biol. Chem.* **267**, 19233-19241 (1992).
2. Powell, W.S., Zhang, Y., and Gravel, S. *Biochemistry* **33**, 3927-3933 (1994).
3. O'Flaherty, J.T., Kuroki, M., Nixon, A.B., et al. *J. Biol. Chem.* **271**, 17821-17828 (1996).
4. O'Flaherty, J.T., Kuroki, M., Nixon, A.B., et al. *J. Immunol.* **157**, 336-342 (1996).
5. Hosoi, T., Koguchi, Y., Sugikawa, E., et al. *J. Biol. Chem.* **277** (35), 31459-31465 (2002).
6. Jones, C.E., Holden, S., Tenailon, L., et al. *Mol. Pharmacol.* **63**(3), 471-477 (2003).

#### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

#### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

#### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 08/30/2022

#### CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA

**PHONE:** [800] 364-9897  
[734] 971-3335

**FAX:** [734] 971-3640

CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM