



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

PRODUCT INFORMATION



Pantoprazole sulfone N-oxide

Item No. 35199

CAS Registry No.: 953787-55-8
Formal Name: 6-(difluoromethoxy)-2-[[[(3,4-dimethoxy-1-oxido-2-pyridinyl)methyl]sulfonyl]-1H-benzimidazole

MF: C₁₆H₁₅F₂N₃O₆S

FW: 415.4


Purity: ≥98%

UV/Vis.: λ_{max}: 279 nm

Supplied as: A crystalline solid

Storage: -20°C

Stability: ≥4 years

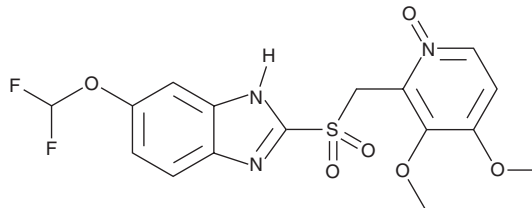

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

Pantoprazole sulfone N-oxide is supplied as a crystalline solid. A stock solution may be made by dissolving the pantoprazole sulfone N-oxide in the solvent of choice, which should be purged with an inert gas. Pantoprazole sulfone N-oxide is soluble in organic solvents such as DMSO and dimethyl formamide. The solubility of pantoprazole sulfone N-oxide in these solvents is approximately 3 and 13 mg/ml, respectively. Pantoprazole sulfone N-oxide is also slightly soluble in ethanol.

Further dilutions of the stock solution into aqueous buffers or isotonic saline should be made prior to performing biological experiments. Ensure that the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. Organic solvent-free aqueous solutions of pantoprazole sulfone N-oxide can be prepared by directly dissolving the crystalline solid in aqueous buffers. The solubility of pantoprazole sulfone N-oxide in PBS (pH 7.2) is approximately 0.5 mg/ml. We do not recommend storing the aqueous solution for more than one day.

Description

Pantoprazole sulfone N-oxide is a potential impurity found in bulk preparations of the proton pump inhibitor pantoprazole.¹

Reference

1. Reddy, G.M., Bhaskar, B.V., Reddy, P.P., *et al.* Structural identification and characterization of potential impurities of pantoprazole sodium. *J. Pharm. Biomed. Anal.* **45(2)**, 201-210 (2007).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 12/12/2022

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM