



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



3,3'-Dipropylthiadicarbocyanine (iodide)

Item No. 36308

CAS Registry No.: 53213-94-8
Formal Name: 3-propyl-2-[5-(3-propyl-2(3H)-benzothiazolylidene)-1,3-pentadien-1-yl]-benzothiazolium, monoiodide
Synonyms: DiS-C₃-(5), DiSC₃(5)
MF: C₂₅H₂₇N₂S₂ • I
FW: 546.5
Purity: ≥98%
Ex./Em. Max: 622/670 nm
UV/Vis.: λ_{max}: 358 nm
Supplied as: A solid
Storage: -20°C
Stability: ≥4 years

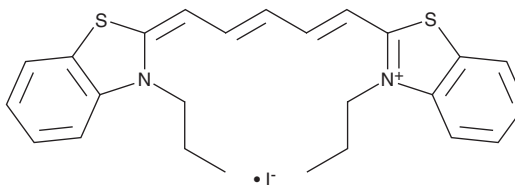

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

3,3'-Dipropylthiadicarbocyanine (iodide) (diSC₃(5) (iodide)) is supplied as a solid. A stock solution may be made by dissolving the diSC₃(5) (iodide) in the solvent of choice, which should be purged with an inert gas. diSC₃(5) (iodide) is soluble in the organic solvent DMSO at a concentration of approximately 2 mg/ml.

Description

diSC₃(5) is a fluorescent dye that has been used to monitor cell membrane potential.¹⁻³ It displays excitation/emission maxima of 622/670 nm, respectively, and, upon cell hyperpolarization, it enters cells, and exhibits a shift in emission maxima to 688 nm and a decrease in fluorescence intensity.¹ When the cell membrane is depolarized, the fluorescence intensity of diSC₃(5) increases as it exits the cells. However, diSC₃(5) (3 μM) also induces membrane hyperpolarization, increases intracellular sodium levels, and decreases intracellular potassium levels in Ehrlich ascites tumor cells, as well as increases potassium membrane permeability by 64.6% in Ehrlich ascites tumor cells in the presence of the Na⁺/K⁺-ATPase inhibitor ouabain.² It also reduces cellular ATP levels in the absence of glucose and induces acid production and inhibits oxygen consumption in the presence of glucose in Ehrlich ascites tumor cells when used at a concentration of 3 μM.³

References

1. Sims, P.J., Waggoner, A.S., Wang, C.H., *et al.* Studies on the mechanism by which cyanine dyes measure membrane potential in red blood cells and phosphatidylcholine vesicles. *Biochemistry* **13**(16), 3315-3330 (1974).
2. Smith, T.C. and Robinson, S.C. The effect of the fluorescent probe, 3,3'-dipropylthiadicarbocyanine iodide, on the membrane potential of Ehrlich ascites tumor cells. *Biochem. Biophys. Res. Commun.* **95**(2), 722-727 (1980).
3. Smith, T.C., Herlihy, J.T., and Robinson, S.C. The effect of the fluorescent probe, 3,3'-dipropylthiadicarbocyanine iodide, on the energy metabolism of Ehrlich ascites tumor cells. *J. Biol. Chem.* **256**(3), 1108-1110 (1981).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 1/21/20221

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM

WWW.CAYMANCHEM.COM