



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



4A3-SC8

Item No. 38155

CAS Registry No.: 1857340-78-3
MF: C₇₅H₁₃₉N₃O₁₆S₄
FW: 1,467.2
Purity: ≥95%
Supplied as: A solution in ethanol
Storage: -20°C
Stability: ≥2 years

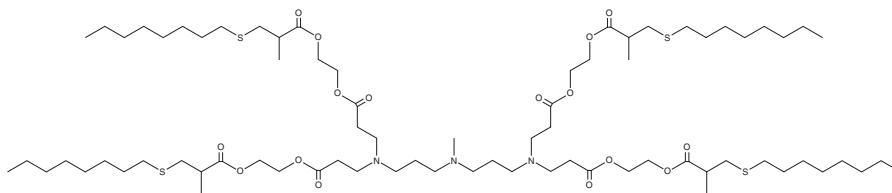

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

4A3-SC8 is supplied as a solution in ethanol. To change the solvent, simply evaporate the ethanol under a gentle stream of nitrogen and immediately add the solvent of choice. 4A3-SC8 is soluble in the organic solvent chloroform.

Description

4A3-SC8 is an ionizable cationic lipid ($pK_a = 6.66$) and a dendrimer that has been used in the generation of lipid nanoparticles (LNPs) for the delivery of siRNA, mRNA, single-guide RNA (sgRNA), ssDNA, and PEGylated BODIPY dyes.¹⁻⁴ LNPs containing 4A3-SC8 and encapsulating siRNA targeting Factor VII selectively accumulate in the liver and decrease hepatic Factor VII activity in mice.³ Intratumoral administration of LNPs containing 4A3-SC8 and encapsulating homologous double-strand repair (HDR) elements correct the DNA encoding GFP with a tyrosine-to-histidine mutation at position 66 (GFP^{Y66H}) in a HEK293 mouse xenograft model.³ Intravenous administration of LNPs containing 4A3-SC8 and encapsulating an mRNA reporter and a PEGylated BODIPY dye have been used for tumor imaging in mice.⁴

References

1. Zhou, K., Nguyen, L.H., Miller, J.B., *et al.* Modular degradable dendrimers enable small RNAs to extend survival in an aggressive liver cancer model. *Proc. Natl. Acad. Sci. USA* **113**(3), 520-525 (2016).
2. Álvarez-Benedicto, E., Farbiak, L., Ramírez, M.M., *et al.* Optimization of phospholipid chemistry for improved lipid nanoparticle (LNP) delivery of messenger RNA (mRNA). *Biomater. Sci.* **10**(2), 549-559 (2022).
3. Farbiak, L., Cheng, Q., Wei, T., *et al.* All-in-one dendrimer-based lipid nanoparticles enable precise HDR-mediated gene editing in vivo. *Adv. Mater.* **33**(30), e2006619 (2021).
4. Xiong, H., Liu, S., Wei, T., *et al.* Theranostic dendrimer-based lipid nanoparticles containing PEGylated BODIPY dyes for tumor imaging and systemic mRNA delivery in vivo. *J. Control. Release* **325**, 198-205 (2020).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 08/20/2024

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM