

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



PRODUCT INFORMATION



C16 Ganglioside G_{M3} - d_9 (d18:1/16:0- d_9) (ammonium salt, synthetic)

Item No. 39225

Formal Name: N-[(1S,2R,3E)-1-[[[O-(N-acetyl- α -neuraminosyl)-(2 \rightarrow 3)-

> O-β-D-galactopyranosyl- $(1\rightarrow 4)$ -β-D-glucopyranosyl]oxy] methyl]-2-hydroxy-3-heptadecen-1-yl]-hexadecanamide-

13,13,14,14,15,15,16,16,16-d₉, monoammonium salt

C16 G_{M3}-d₉, N-Hexadecanoyl-d₉ Synonyms:

(13,13,14,14,15,15,16,16,16)-Monosialoganglioside G_{M3}

MF: C₅₇H₉₄D₉N₂O₂₁ • NH₄

FW: 1,179.5


Chemical Purity: ≥95% (C16 Ganglioside G_{M3} (d18:1/16:0) (ammonium

Deuterium

Incorporation: ≥99% deuterated forms (d₁-d₉); ≤1% d₀

Supplied as: A film Storage: -20°C Stability: ≥2 years

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

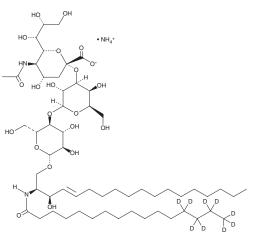

Laboratory Procedures

C16 Ganglioside G_{M3} - d_9 (d18:1/16:0- d_9) is intended for use as an internal standard for the quantification of C16 ganglioside G_{M3} by GC- or LC-MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

C16 Ganglioside G_{M3} -d₉ (d18:1/16:0-d₉) is supplied as a solid. A stock solution may be made by dissolving the C16 ganglioside G_{M3} - d_9 (d18:1/16:0- d_9) in the solvent of choice, which should be purged with an inert gas. C16 Ganglioside G_{M3}-d₉ (d18:1/16:0-d₉) is soluble in a 2:1:0.1 solution of chloroform:methanol:DI water.

Description

C16 Ganglioside G_{M3} -d $_9$ (d18:1/16:0-d $_9$) is intended for use as an internal standard for the quantification of C16 ganglioside G_{M3} by GC- or LC-MS. C16 Ganglioside G_{M3} (d18:1/16:0) is a monosialylated ganglioside. It has been found in human fetal frontal lobe tissue isolated in the 27th gestational week but not in human fetal anencephalic residual brain tissue isolated in the 28^{th} gestational week. C16 Ganglioside G_{M3} (d18:1/16:0) is higher in the serum of patients with visceral fat accumulation (VFA) and patients with VFA and dyslipidemia compared to healthy lean control subjects. 1 Serum levels of C16 ganglioside G_{M3} (d18:1/16:0) are inversely correlated with erythrocyte counts in patients with hematological diseases, including patients with lymphoid neoplasms.3

References

- 1. Almeida, R., Mosoarca, C., Chirita, M., et al. Anal. Biochem. 378(1), 43-52 (2008).
- 2. Veillon, L., Go, S., Matsuyama, W., et al. PLoS One 10(6), e0129645 (2015).
- Nishikawa, M., Kurano, M., Nitta, T., et al. Sci. Rep. 9(1), 6308 (2019).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website

Copyright Cayman Chemical Company, 09/08/2023

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD ANN ARBOR, MI 48108 · USA PHONE: [800] 364-9897

[734] 971-3335

FAX: [734] 971-3640 CUSTSERV@CAYMANCHEM.COM WWW.**CAYMANCHEM**.COM