



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)

# PRODUCT INFORMATION



## LipidLaunch™ LNP SM-102 (Luciferase)


Item No. 39318

### Overview and Properties

**Storage:** -80°C (as supplied)  
**Stability:** ≥6 months  
**Supplied in:** PBS, pH 7.4, with 10% sucrose

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

### Images

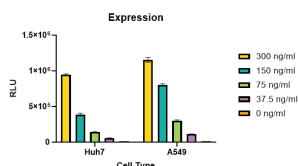

Transfection of cells with luciferase-encapsulating SM-102 particles, Huh7 hepatocytes (left), and A549 lung epithelial cells (right), which were incubated with Cayman's LipidLaunch™ LNP SM-102 (Luciferase) (Item No. 39318) at indicated RNA concentrations for 24 hours. Luciferase expression (RLU) is shown.

Parameter	Range/Value
Size	75-125 nm
Polydispersity index (PDI)	<0.2
Encapsulation efficiency (%EE)	>95%
mRNA concentration	Batch specific, 25-100 µg/ml
mRNA/vial	5 µg

### Description

LipidLaunch™ LNP SM-102 (Luciferase) is a solution containing lipid nanoparticles (LNPs) composed of the ionizable cationic amino lipid SM-102 (Item No. 33474), cholesterol (Item No. 9003100), the phospholipid 1,2-distearoyl-*sn*-glycero-3-PC (Item Nos. 15100 | 39189), and the lipid excipient DMG-PEG(2000) (Item No. 33945) at a molar ratio of 50:38.5:10:1.5 and firefly luciferase protein. It is intended for proof-of-concept experiments to determine whether SM-102-based LNPs can effectively lead to expression of a protein of interest in a target cell type, either *in vitro* or *in vivo*. Firefly luciferase is an enzyme that catalyzes the oxidation of ATP-dependent D-luciferin, resulting in chemiluminescence at a wavelength of approximately 560 nm.

**Suggested *in vitro* use:** Thaw LNPs on ice with occasional gentle swirling (*do not vortex*). Using a gentle pipetting technique, dilute 1:100-1:500 in complete cell culture media (with serum) and add to subconfluent cells in a luminescence-compatible tissue culture plate. Use a luciferase reporter assay substrate and read luminescence on a plate reader. Optimal conditions are highly dependent on cell type.

**WARNING**  
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

**SAFETY DATA**  
This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

**WARRANTY AND LIMITATION OF REMEDY**  
Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 09/14/2023

**CAYMAN CHEMICAL**  
1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA  
PHONE: [800] 364-9897  
[734] 971-3335  
FAX: [734] 971-3640  
CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM