



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


PRODUCT INFORMATION



C18

Item No. 39432

CAS Registry No.: 2922824-09-5
Formal Name: 2-chloro-N-(2,3-dihydro-1,4-benzodioxin-6-yl)-N-(4-fluoro-3,5-dimethoxyphenyl)-acetamide
MF: C₁₈H₁₇ClFNO₅
FW: 381.8
Purity: ≥98%
Supplied as: A solid
Storage: -20°C
Stability: ≥3 years

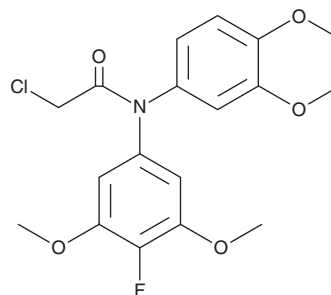

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

C18 is supplied as a solid. A stock solution may be made by dissolving the C18 in the solvent of choice, which should be purged with an inert gas. C18 is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide (DMF). The solubility of C18 in DMSO and DMF is approximately 1 mg/ml. C18 is slightly soluble in ethanol.

Description

C18 is a covalent inhibitor of glutathione peroxidase 4 (GPX4; IC₅₀ = 120 nM).¹ It inhibits the growth of MDA-MB-468, BT-549, and MDA-MB-231 breast cancer cells (IC₅₀s = 10, 75, and 12 nM, respectively). C18 (10 and 20 μM) increases lipid peroxidation and oxidized malondialdehyde (MDA) levels in MDA-MB-231 cells, an effect that can be prevented by the ferroptosis inhibitor ferrostatin-1 (Item No. 17729). It decreases mitochondrial volume and induces degradation of mitochondrial crests and outer membranes in MDA-MB-231 cells when used at a concentration of 20 μM. *In vivo*, C18 (5, 10, or 20 mg/kg) reduces tumor volume and weight, without affecting body weight, in an MDA-MB-231 mouse xenograft model.

Reference

1. Chen, T., Leng, J., Tan, J., *et al.* Discovery of novel potent covalent glutathione peroxidase 4 inhibitors as highly selective ferroptosis inducers for the treatment of triple-negative breast cancer. *J. Med. Chem.* **66**(14), 10036-10059 (2023).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 09/14/2023

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM