



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)


[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)

# PRODUCT INFORMATION



**XIE62-1004**  
Item No. 40446

**CAS Registry No.:** 2421146-32-7  
**Formal Name:** 2-[[[3,4-bis(phenylmethoxy)phenyl]methyl]amino]-ethanol, monohydrochloride  
**MF:** C<sub>23</sub>H<sub>25</sub>NO<sub>3</sub> • HCl  
**FW:** 399.9  
**Purity:** ≥98%  
**Supplied as:** A solid  
**Storage:** -20°C  
**Stability:** ≥4 years

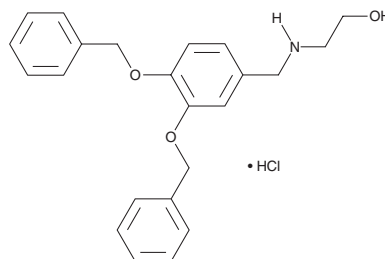

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

## Laboratory Procedures

XIE62-1004 is supplied as a solid. A stock solution may be made by dissolving the XIE62-1004 in the solvent of choice, which should be purged with an inert gas. XIE62-1004 is soluble in DMSO. XIE62-1004 is slightly soluble in acetonitrile.

XIE62-1004 is slightly soluble in aqueous solutions. To enhance aqueous solubility, dilute the organic solvent solution into aqueous buffers or isotonic saline. If performing biological experiments, ensure the residual amount of organic solvent is insignificant, since organic solvents may have physiological effects at low concentrations. We do not recommend storing the aqueous solution for more than one day.

## Description

XIE62-1004 is an inducer of autophagy.<sup>1</sup> It increases the number of p62 aggregates, as well as the number of autophagosomes and autolysosomes, in wild-type HeLa cervical cancer cells, but not in HeLa cells deficient in the gene encoding p62, *SQSTM1*, when used at a concentration of 10 μM. XIE62-1004 induces autophagic flux and the conversion of LC3 autophagy marker I (LC3-I) to LC3-II in HeLa cells. It increases the degradation of mutant huntingtin (mHTT<sup>Q103</sup>) protein expressed in HeLa cells.

## Reference

1. Cha-Molstad, H., Yu, J.E., Feng, Z., *et al.* p62/SQSTM1/Sequestosome-1 is an N-recogin of the N-end rule pathway which modulates autophagosome biogenesis. *Nat. Commun.* **8(1)**, 102 (2017).

### WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

### SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

### WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 03/05/2024

## CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD  
ANN ARBOR, MI 48108 · USA

**PHONE:** [800] 364-9897  
[734] 971-3335

**FAX:** [734] 971-3640

CUSTSERV@CAYMANCHEM.COM  
WWW.CAYMANCHEM.COM