



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



## 4E-BP1 (Phospho-Thr<sup>37/46</sup>) TR-FRET Assay Kit

*Powered by Bioauxilium's THUNDER™ TR-FRET Technology*

---

Item No. 500224

[www.caymanchem.com](http://www.caymanchem.com)

Customer Service 800.364.9897

Technical Support 888.526.5351

1180 E. Ellsworth Rd • Ann Arbor, MI • USA

---

## TABLE OF CONTENTS

GENERAL INFORMATION	3	Materials Supplied
	3	Safety Data
	4	Precautions
	4	Before You Start
	5	If You Have Problems
	5	Storage and Stability
	6	THUNDER™ General Information
INTRODUCTION	7	Materials Needed but Not Supplied
	8	Background
	9	About This Assay
PRE-ASSAY PREPARATION	10	Principle Of This Assay
	11	Assay Optimization
	12	Reagent Preparation
ASSAY PROTOCOL	15	TR-FRET Plate Reader Settings
	16	Workflow
	17	Assay Summary
ANALYSIS	18	Performing the Assay: 2-plate (Transfer) Protocol
	25	Performing the Assay: 1-plate (All in One Well) Protocol
	30	Calculations
RESOURCES	31	Performance Characteristics
	35	Troubleshooting
	37	References
	38	Notes
	39	Warranty and Limitation of Remedy

## GENERAL INFORMATION

### Materials Supplied

Item Number	Item Name	96 wells Quantity/Size	480 wells Quantity/Size
400271	Europium-Labeled 4E-BP1 (Phospho-Thr <sup>37/46</sup> ) Antibody	1 vial/5 µl	1 vial/25 µl
400272	Acceptor-Labeled 4E-BP1 (Phospho-Thr <sup>37/46</sup> ) Antibody	1 vial/20 µl	1 vial/100 µl
400211	Lysis Buffer 2 (5X)	1 vial/1 ml	4 vials/1.25 ml
400225	Detection Buffer (10X)	1 vial/50 µl	1 vial/250 µl
400297	4E-BP1 Positive Control Lysate	1 vial/100 µl	1 vial/200 µl
400257	Phosphatase Inhibitor Cocktail (100X)	1 vial/50 µl	1 vial/250 µl

If any of the items listed above are damaged or missing, please contact our Customer Service department at (800) 364-9897 or (734) 971-3335. We cannot accept any returns without prior authorization.



**WARNING: THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.**

### Safety Data

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent *via* email to your institution.

## Precautions

**Please read these instructions carefully before beginning this assay.**

Do not mix or substitute reagents or materials from other kit lots or kits. Kits are quality control tested as a set of components and performance cannot be guaranteed if utilized separately or substituted.

We cannot guarantee the performance of the product outside the conditions detailed in this kit booklet.

The kits are designed for the detection of endogenous cellular proteins across a wide variety of cell lines. However, until each cell line in particular is tested, the possibility of the presence of undetectable levels of the target protein cannot be excluded.

Users should ensure that their cell line has measurable levels of the target protein. Expression levels of signaling proteins in different cell types vary widely. The cell line used for the assay validation of this kit is shown in the figures starting on page 31.

## Before You Start

**Please note the following:**

ONLY white plates should be used for TR-FRET.

DO NOT modify the assay protocol or volumes.

DO optimize the cell density, serum starvation (optional), and stimulation or inhibition parameters.

ALWAYS use the included positive control lysate for every assay.

## If You Have Problems

### Technical Service Contact Information

Phone: 888-526-5351 (USA and Canada only) or 734-975-3888

Email: techserv@caymanchem.com

In order for our staff to assist you quickly and efficiently, please be ready to supply the lot number of the kit (found on the outside of the box).

## Storage and Stability

This kit will perform as specified if stored as directed at -80°C and used before the expiration date indicated on the outside of the box.

## THUNDER™ General Information

THUNDER™ TR-FRET Cell Signaling Assay Kits are designed for the semi-quantitative measurement of phosphorylated and/or total (both phosphorylated and unphosphorylated) proteins in cell lysates using homogeneous (no wash) TR-FRET technology. The kits are compatible with both adherent and suspended cells.

THUNDER™ TR-FRET Cell Signaling Assay Kits are based on Bioauxilium's enhanced proprietary time-resolved Förster resonance energy transfer (TR-FRET) technology. THUNDER™ assays can be read on most commercially available TR-FRET-compatible plate readers (a list of suitable TR-FRET readers can be found at [www.Bioauxilium.com](http://www.Bioauxilium.com)). TR-FRET-based assays are homogeneous because they do not require any washing or separation steps. In addition, the THUNDER™ assays use a standardized, simple, and rapid "add-incubate-measure" protocol with a single step reagent addition. This streamlined assay protocol dramatically decreases hands-on time and provides a powerful alternative to cumbersome, error-prone and time-consuming techniques such as Western blot and ELISA.

THUNDER™ TR-FRET Cell Signaling Assay Kits contain the essential reagents necessary to carry out the measurement of signaling proteins in cells, with the exception of the plate(s).

## Materials Needed But Not Supplied

1. A plate reader equipped with a TR-FRET option
2. Adjustable pipettes; multichannel or repeating pipettor recommended
3. A source of ultrapure water, with a resistivity of 18.2 MΩ-cm and total organic carbon (TOC) levels of <10 ppb, is recommended. Pure water - glass-distilled or deionized - may not be acceptable. *NOTE: UltraPure Water is available for purchase from Cayman (Item No. 400000).*
4. Culture plate: 96-well clear, flat-bottom polystyrene tissue culture-treated plate(s) for culturing cells when using the 2-plate (transfer) assay protocol. *NOTE: Do not use this type of plate for the 1-plate (all in one well) assay protocol.*
5. Detection plate (96-well plate option): Half-area, 96-well white plate(s) for TR-FRET detection when using the 2-plate (transfer) assay protocol
6. Detection plate (384-well plate option): Low-volume, 384-well white plate(s) for TR-FRET detection when using the 1- or 2-plate assay protocols
7. Adhesive sealing film for plates
8. Orbital microplate shaker

## INTRODUCTION

### Background


Eukaryotic translation initiation factor 4E (eIF4E)-binding protein 1 (4E-BP1) is a member of the eIF4E-binding protein family and a repressor of translation.<sup>1</sup> Non-phosphorylated 4E-BP1 inhibits the initiation of mRNA translation by binding to eIF4E and preventing it from being incorporated into the 48S pre-initiation complex required for the recruitment of 40S ribosomal subunits to the 5' end of mRNA.<sup>1,2</sup> It is inactivated by phosphorylation at positions Thr<sup>37</sup> and Thr<sup>46</sup> (phospho-Thr<sup>37/46</sup>) by mammalian target of rapamycin complex 1 (mTORC1), or in an mTOR-independent manner by glycogen synthase kinase 3 $\beta$  (GSK3 $\beta$ ) or leucine-rich repeat kinase 2 (LRRK2), and detaches from eIF4E, allowing translation to proceed.<sup>2</sup> 4E-BP1 has a tumor suppressing role in cancer cells by reducing the overall rate of translation and inhibiting the translation of certain transcripts, such as cell cycle progression and angiogenic factors.<sup>1</sup> However, its role can become pro-oncogenic under conditions of cellular stress in the tumor microenvironment by inducing a switch to cap-independent translation and increasing the translation of factors that are proangiogenic under hypoxic conditions. Cytoplasmic levels of phosphorylated (inactive) 4E-BP1 are negatively correlated with recurrence-free survival in patients with progesterone receptor-expressing breast cancer, while nuclear levels of non-phosphorylated (active) 4E-BP1 are positively correlated with recurrence-free survival.<sup>3</sup> High levels of 4E-BP1 (phospho-Thr<sup>37/46</sup>) are associated with poor prognosis in patients with early stage esophageal squamous cell carcinoma.<sup>4</sup> 4E-BP1 (phospho-Thr<sup>37/46</sup>) levels have also been used as a marker of protein synthesis.<sup>5,6</sup>

### About This Assay

This 4E-BP1 (Phospho-Thr<sup>37/46</sup>) TR-FRET Assay Kit uses a homogeneous TR-FRET assay method amenable to rapid measurement of phosphorylated protein levels in cells. This 4E-BP1 (Phospho-Thr<sup>37/46</sup>) TR-FRET Assay Kit is robust, with a Z' value of 0.73, and suitable for screening a large number of samples. The signal is stable at room temperature for at least 24 hours, affording flexibility in read times. The amount of reagents provided is sufficient for testing either 96 or 480 phosphorylated protein wells, depending on the size of the kit.

## Principle Of This Assay

This assay is based on the traditional sandwich immunoassay principle (Figure 1, below). Following cell treatment, cells are lysed with the specific lysis buffer provided in the kit. Then, the target protein in the cell lysates is detected with a pair of fluorophore-labeled antibodies reactive to human samples.

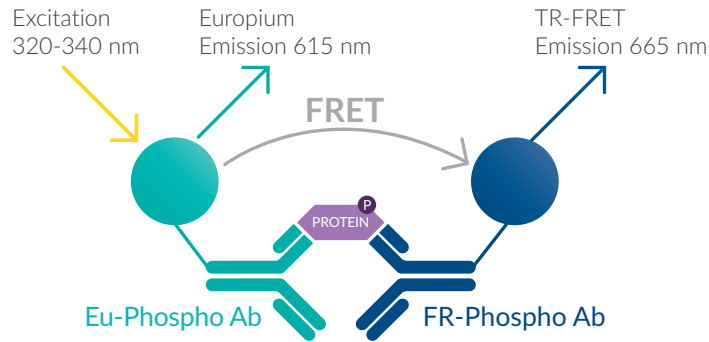

**Figure 1. Schematic of the TR-FRET cell signaling assay principle**

The first antibody is labeled with a long-lifetime donor fluorophore (a europium chelate; Eu-Phospho Ab) and the second with a far-red acceptor fluorophore (FR-Phospho Ab). The binding of the two labeled antibodies to distinct epitopes on the target protein takes place in solution and brings the two dyes into close proximity. Excitation of the donor Eu chelate molecules with a flash lamp (320 or 340 nm) or laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which, in turn, emit a TR-FRET signal at 665 nm. The signal at 665 nm is proportional to the concentration of target protein in the cell lysate. Residual energy from the Eu chelate generates light at 615 nm, which can be used as an internal standard to normalize light emitted at 665 nm.

TR-FRET assays exhibit very low background fluorescence levels and high signal-to-background (S/B) ratios. The data can be expressed and analyzed as either the signal at 665 nm or the 665 nm/615 nm ratio. The ratiometric measurement further increases assay reproducibility and robustness.

## PRE-ASSAY PREPARATION

### Assay Optimization

A critical step in performing any cell-based assay is the optimization of cell culture and treatment conditions. The following protocol assumes that both the cell number and treatment conditions have been previously optimized, as these key parameters often vary for each cell line. It is, therefore, strongly recommended to optimize these parameters in order to maximize the assay signal and ensure optimum performance with a high S/B ratio.

Cell number, serum starvation (optional), and stimulation or inhibition time (at either room temperature or 37°C) should be optimized for each cell line and target protein. Cell numbers that are too high or too low can negatively influence the activation of intracellular signaling pathways. Cell seeding densities of 40,000-80,000 cells/well for adherent cells or 100,000-200,000 cells/well for suspended cells are generally acceptable for most cell lines. Of note, the optimal length of time for stimulation can vary widely among cell lines from a few minutes to more than one hour. As such, a time-course study is strongly recommended to determine the optimal stimulation time, ideally at both room temperature and 37°C, since incubation temperature has an effect on the kinetics of target protein stimulation. Additional assay development guidelines are available on Bioauxilium's website ([www.Bioauxilium.com](http://www.Bioauxilium.com)).

## Reagent Preparation

The instructions described below are for testing the entire number of wells in each kit. Adjust volumes accordingly when testing fewer wells.

Bring all reagents to room temperature prior to use.

Centrifuge all tubes before use to improve recovery of content (2,000 x g, 10-15 seconds).

Mix the lysis and detection buffers and the positive control lysate by vortexing gently before use. Do NOT vortex the antibodies.

Use ultrapure water (18 M $\Omega$ -cm) to dilute the lysis and detection buffers.

NOTE: It is recommended to test all samples and controls at least in duplicate.

NOTE: ALWAYS include a positive control using the positive control lysate provided.

### 1. Supplemented Lysis Buffer

**Supplemented Lysis Buffer 2 (1X) for the 2-Plate (Transfer) Assay Protocol with Adherent Cells:** The supplemented Lysis Buffer 2 (1X) is designed for use in the 2-plate (transfer) assay protocol using adherent cells (see page 18). Each well requires 50  $\mu$ l of supplemented Lysis Buffer 2 (1X). Dilute the Lysis Buffer 2 (5X) (Item No. 400211) with ultrapure water and add the Phosphatase Inhibitor Cocktail (100X) (Item No. 400257), which contains sodium fluoride (NaF), sodium orthovanadate (Na<sub>3</sub>VO<sub>4</sub>), and glycerophosphate at 100, 200, and 200 mM, respectively, to final NaF, Na<sub>3</sub>VO<sub>4</sub>, and glycerophosphate concentrations of 1, 2, and 2 mM, respectively. NOTE: It is mandatory to supplement Lysis Buffer 2 with the Phosphatase Inhibitor Cocktail (100X). Store unused Lysis Buffer 2 (1X) at 4°C; it will be stable for approximately two days.

OR

**Supplemented Lysis Buffer 2 (5X) for the 2-Plate (Transfer) Assay Protocol with Suspension Cells or the 1-Plate Assay Protocol:** The supplemented Lysis Buffer 2 (5X) is designed for use in the 2-plate (transfer) assay protocol using suspension cells (see page 20) or the 1-plate (all in one well) assay protocol for adherent or suspension cells (see page 25). Each well requires 10  $\mu$ l (transfer) or 3  $\mu$ l (all in one well) of supplemented Lysis Buffer 2 (5X). Directly add the Phosphatase Inhibitor Cocktail (100X) (Item No. 400257), which contains sodium fluoride (NaF), sodium orthovanadate (Na<sub>3</sub>VO<sub>4</sub>), and glycerophosphate at 100, 200, and 200 mM, respectively, to final NaF, Na<sub>3</sub>VO<sub>4</sub>, and glycerophosphate concentrations of 5, 10, and 10 mM, respectively. NOTE: It is mandatory to supplement Lysis Buffer 2 with the Phosphatase Inhibitor Cocktail (100X). Store unused Lysis Buffer 2 (5X) at 4°C; it will be stable for approximately two days.

### 2. Detection Buffer (1X)

Dilute 50  $\mu$ l (96 wells) or 250  $\mu$ l (480 wells) Detection Buffer (10X) (Item No. 400225) with 0.45 ml (96 wells) or 2.25 ml (480 wells), respectively, of ultrapure water. Store unused Detection Buffer (1X) at 4°C; it will be stable for approximately two days.



### 3. Antibody Detection Mix (4X)

Prepare and mix just before use. *NOTE: Due to the low reagent volumes, the antibodies are diluted with Detection Buffer (1X) directly in the vial when assaying only 96 wells.*

*NOTE: Each well requires 5 µl of Antibody Detection Mix (4X).*

**Antibody Detection Mix (4X) (96 wells):** Add 255 µl of Detection Buffer (1X) into the vial containing 5 µl of Europium-Labeled 4E-BP1 (Phospho-Thr<sup>37/46</sup>) Antibody (Eu-Phospho Ab) (Item No. 400271). Add 240 µl of Detection Buffer (1X) into the vial containing 20 µl of Acceptor-Labeled 4E-BP1 (Phospho-Thr<sup>37/46</sup>) Antibody (FR-Phospho Ab) (Item No. 400272). Gently mix the diluted Eu-Phospho Ab and FR-Phospho Ab solutions together.

OR

**Antibody Detection Mix (4X) (480 wells):** Gently mix 25 µl of Eu-Phospho Ab (Item No. 400271) with 1,275 µl of Detection Buffer (1X). Gently mix 100 µl of FR-Phospho Ab (Item No. 400272) with 1,200 µl of Detection Buffer (1X). Gently mix the diluted Eu-Phospho Ab and FR-Phospho Ab solutions together.

Store unused Antibody Detection Mix (4X) working solution at 4°C; it will be stable for approximately two days.

### 4. Positive Control Lysate

The 4E-BP1 Positive Control Lysate (Item No. 400297) is supplied ready to use. The thawed positive control lysate can be aliquoted, refrozen at -80°C, and thawed at least three more times.

## TR-FRET Plate Reader Settings

We recommend reading the TR-FRET assays at two wavelengths, detecting both the emission from the Eu chelate donor fluorophore at 615 nm and the acceptor fluorophore at 665 nm. Table 1, below, provides instrument settings to be used as guidelines.

Parameter	TR-FRET-compatible Plate Reader	
	Flash lamp excitation	Laser excitation
Excitation filter	320 nm (or 340 nm)	N/A
Emission filter	615 nm (or 620 nm)	615 nm (or 620 nm)
Delay time	90 µs	50 µs
Flash energy level	100% or High	100%
Number of flashes	100	20
Window (integration time)	300 µs	100 µs

Table 1. Recommended TR-FRET plate reader settings

## ASSAY PROTOCOL

### Workflow

The THUNDER™ TR-FRET Cell Signaling Assay workflow consists of 3 simple steps.


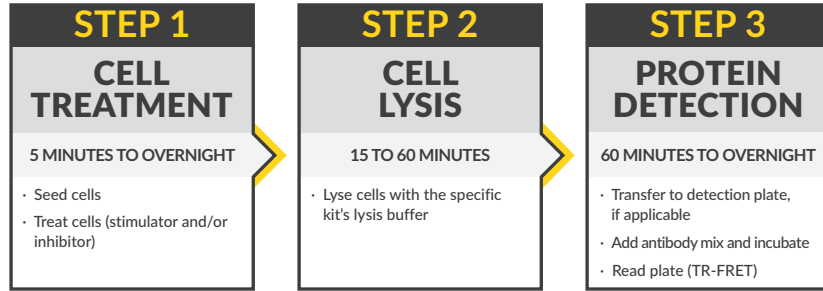


Figure 2. Assay workflow


### Assay Summary

The THUNDER™ TR-FRET Cell Signaling Assays can be run using one of two possible protocols.

2-Plate (transfer) Protocol

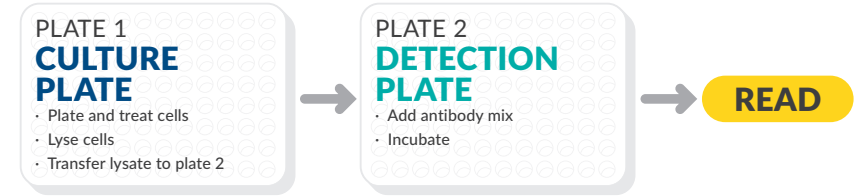

1-Plate (all in one well) Protocol

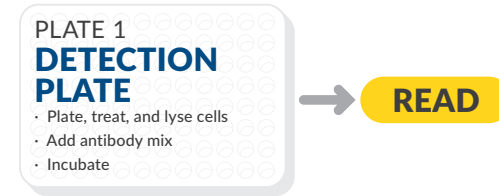

**Figure 3. Assay summary** Each assay uses the same total volume of 20 µl for TR-FRET detection. The 2-plate protocol is conducted in either a half-area 96-well or low-volume 384-well plate(s), while the 1-plate protocol is conducted in a low-volume 384-well plate(s).

## Performing the Assay: 2-plate (Transfer) Protocol

### Adherent Cells

This is a transfer protocol that is conducted in two different types of plates: cell culture and lysis are conducted in a 96-well culture plate(s), whereas detection is conducted in either a white, half-area 96-well assay plate(s) or a white, low-volume 384-well assay plate(s), with a volume of 20  $\mu\text{l}$  per well for TR-FRET detection. A summary of this protocol is provided in Table 2 (see page 22).

### Cell Treatment

1. Dispense 50  $\mu\text{l}$  of cells at the pre-optimized density into a 96-well tissue culture-treated plate(s) in appropriate culture medium.
2. Incubate overnight at 37°C and 5% CO<sub>2</sub>.
3. Cell Stimulation or Inhibition

**Stimulation:** Add 50  $\mu\text{l}$  of stimulator (2X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperature needs to be determined.*

**No treatment:** Add 50  $\mu\text{l}$  of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells.

OR

**Inhibition:** Add 25  $\mu\text{l}$  of inhibitor (4X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. Add 25  $\mu\text{l}$  of stimulator (4X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperatures need to be determined.*

**No treatment:** Add 25  $\mu\text{l}$  of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells. Add another 25  $\mu\text{l}$  of serum-free medium and incubate for the same amount of time, and at the same temperature, as treated cells.

### Cell Lysis

1. Carefully remove the cell culture medium by aspirating the supernatant.
2. Immediately add 50  $\mu\text{l}$  of supplemented Lysis Buffer 2 (1X). The volume of supplemented Lysis Buffer 2 (1X) may be optimized from 25-100  $\mu\text{l}$ .
3. Incubate for 30 minutes at room temperature on an orbital plate shaker at 400 rpm. The incubation time in supplemented Lysis Buffer 2 (1X) may be optimized from 15-60 minutes.

*NOTE: If samples cannot be measured immediately, store at -80°C.*

### TR-FRET Detection

1. Carefully pipette 15  $\mu\text{l}$  of cell lysate from the 96-well culture plate to a well of either a white, half-area 96-well or a white, low-volume 384-well plate.
2. It is recommended to add 15  $\mu\text{l}$  of 4E-BP1 Positive Control Lysate (undiluted) and 15  $\mu\text{l}$  of Lysis Buffer 2 (1X) (negative control) to separate assay wells.
3. Add 5  $\mu\text{l}$  of Antibody Detection Mix (4X) to each of the assay wells.
4. Cover the plate(s) with a plate sealer and incubate for 4 hours at room temperature.
5. Gently remove the adhesive plate sealer. Read on a TR-FRET-compatible plate reader. *NOTE: The same plate can be read several times without a negative effect on the assay performance, and the signal is stable for at least 24 hours at room temperature.*

*NOTE: A summary of this pipetting protocol is provided in Table 4 (see page 24).*

## Suspension Cells

This is a transfer protocol that is conducted in two different types of plates: cell culture and lysis are conducted in a 96-well culture plate(s), whereas detection is conducted in either a white, half-area 96-well assay plate(s) or a white, low-volume 384-well assay plate(s), with a volume of 20  $\mu$ l per well for TR-FRET detection. A summary of this protocol is provided in Table 3 (see page 23).

### Cell Treatment

1. Dispense 20  $\mu$ l of cells at the pre-optimized density into a 96-well tissue culture-treated plate(s) in appropriate culture medium.
2. Directly proceed to cell stimulation or inhibition or incubate 2-4 hours at 37°C and 5% CO<sub>2</sub> prior to stimulation or inhibition. This step may be optimized.
3. Cell Stimulation or Inhibition

**Stimulation:** Add 20  $\mu$ l of stimulator (2X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperature needs to be determined.*

**No treatment:** Add 20  $\mu$ l of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells.

OR

**Inhibition:** Add 10  $\mu$ l of inhibitor (4X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. Add 10  $\mu$ l of stimulator (4X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperatures need to be determined.*

**No treatment:** Add 10  $\mu$ l of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells. Add another 10  $\mu$ l of serum-free medium and incubate for the same amount of time, and at the same temperature, as treated cells.

## Cell Lysis

1. Add 10  $\mu$ l supplemented Lysis Buffer 2 (5X).
2. Incubate for 30 minutes at room temperature on an orbital plate shaker at 400 rpm. The incubation time in supplemented Lysis Buffer 2 (5X) may be optimized from 15-60 minutes.

*NOTE: If samples cannot be measured immediately, store at -80°C.*

### TR-FRET Detection

Following cell lysis, proceed to the TR-FRET detection step as described for the standard 2-plate (transfer) protocol for adherent cells (see page 19).

2-Plate (Transfer) Assay Summary				
Step	Adherent Cells			
Cell Treatment	Stimulation	No Treatment	Inhibition	No Treatment
	50 µl cells	50 µl cells	50 µl cells	50 µl cells
	Incubate cells overnight			
	50 µl stimulator (2X)	50 µl serum-free medium	25 µl inhibitor (4X)	25 µl serum-free medium
	Incubate for pre-optimized time			
			25 µl stimulator (4X)	25 µl serum-free medium
			Incubate for pre-optimized time	
Cell Lysis	Remove media			
	50 µl supplemented Lysis Buffer 2 (1X)*			
	Incubate 30 minutes on an orbital shaker			
Protein Detection	15 µl lysate			
	5 µl Antibody Detection Mix (4X)			
	Cover and incubate 4 hours			
	Read TR-FRET signal			

**Table 2. Assay summary for the 2-plate (transfer) protocol with adherent cells**  
\*The lysis buffer must be supplemented with the Phosphatase Inhibitor Cocktail (100X) as described in the Reagent Preparation section.

2-Plate (Transfer) Assay Summary				
Step	Suspension Cells			
Cell Treatment	Stimulation	No Treatment	Inhibition	No Treatment
	20 µl cells	20 µl cells	20 µl cells	20 µl cells
	20 µl stimulator (2X) OR incubate for 2-4 hours, then add 20 µl stimulator (2X)	20 µl serum-free medium OR incubate for 2-4 hours, then add 20 µl serum-free medium	10 µl inhibitor (4X) OR incubate for 2-4 hours, then add 10 µl inhibitor (4X)	10 µl serum-free medium OR incubate for 2-4 hours, then add 10 µl serum-free medium
	Incubate for pre-optimized time			
			10 µl inhibitor (4X)	10 µl serum-free medium
			Incubate for pre-optimized time	
Cell Lysis	10 µl supplemented Lysis Buffer 2 (5X)*			
	Incubate 30 minutes on an orbital shaker			
Protein Detection	15 µl lysate			
	5 µl Antibody Detection Mix (4X)			
	Cover and incubate 4 hours			
	Read TR-FRET signal			

**Table 3. Assay summary for the 2-plate (transfer) protocol with suspension cells**  
\*The lysis buffer must be supplemented with the Phosphatase Inhibitor Cocktail (100X) as described in the Reagent Preparation section.

	Untreated Cells	Treated Cells	Positive Control Lysate	Negative Control Lysate
Cell lysate (untreated cells)	15 $\mu$ l			
Cell lysate (treated cells)		15 $\mu$ l		
Positive Control Lysate			15 $\mu$ l	
Lysis Buffer 2 (1X)				15 $\mu$ l
Antibody Detection Mix (4X)	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l
<b>Total assay volume</b>	<b>20 <math>\mu</math>l</b>	<b>20 <math>\mu</math>l</b>	<b>20 <math>\mu</math>l</b>	<b>20 <math>\mu</math>l</b>

**Table 4. Summary of pipetting protocol 2-Plate (transfer) protocol after lysis and prior to TR-FRET**

## Performing the Assay: 1-plate (All in One Well) Protocol

### Adherent and Suspension Cells

This is an all-in-one-well protocol. No transfer step is needed. Conduct the assay in a white, low-volume 384-well assay plate(s) with a total assay volume of 20  $\mu$ l per well. A summary of this protocol is provided in Tables 5-7 (see pages 27-29).

#### Cell Treatment

1. Dispense 8  $\mu$ l of cells at the pre-optimized density in serum-free medium into a white, low-volume 384-well assay plate(s).
2. Cell Stimulation or Inhibition

**Stimulation:** Add 4  $\mu$ l of stimulator (3X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperature needs to be determined.*

**No treatment:** Add 4  $\mu$ l of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells.

**OR**

**Inhibition:** Add 2  $\mu$ l of inhibitor (6X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. Add 2  $\mu$ l of stimulator (6X) diluted in serum-free medium. Incubate for the pre-optimized time at either room temperature or 37°C. *NOTE: Optimal incubation temperatures need to be determined.*

**No treatment:** Add 2  $\mu$ l of serum-free medium. Incubate for the same amount of time, and at the same temperature, as treated cells. Add another 2  $\mu$ l of serum-free medium and incubate for the same amount of time, and at the same temperature, as treated cells.

## Cell Lysis

1. Add 3  $\mu$ l of supplemented Lysis Buffer 2 (5X).
2. Incubate for 30 minutes at room temperature on an orbital plate shaker at 400 rpm. The incubation time in supplemented Lysis Buffer 2 (5X) may be optimized from 15-60 minutes.

*NOTE: If samples cannot be measured immediately, store at -80°C.*

## TR-FRET Detection

1. Add 15  $\mu$ l of 4E-BP1 Positive Control Lysate (undiluted) and 15  $\mu$ l of Lysis Buffer 2 (1X) (negative control) to separate assay wells.
2. Add 5  $\mu$ l of Antibody Detection Mix (4X) prepared in Detection Buffer (1X) to each of the assay wells.
3. Cover the plate(s) with a plate sealer and incubate for 4 hours at room temperature.
4. Gently remove the adhesive plate sealer. Read on a TR-FRET-compatible plate reader. *NOTE: The same plate can be read several times without a negative effect on the assay performance, and the signal is stable for at least 24 hours at room temperature.*

1-Plate Assay Summary				
Step	Adherent or suspension cells			
White, low-volume 384-well plate(s)				
Cell Treatment	Stimulation	No treatment	Inhibition	No treatment
	8 $\mu$ l cells	8 $\mu$ l cells	8 $\mu$ l cells	8 $\mu$ l cells
	4 $\mu$ l stimulator (3X)	4 $\mu$ l serum-free medium	2 $\mu$ l inhibitor (6X)	2 $\mu$ l serum-free medium
	Incubate for pre-optimized time			
			2 $\mu$ l stimulator (6X)	2 $\mu$ l serum-free medium
			Incubate for pre-optimized time	
Cell Lysis	3 $\mu$ l supplemented Lysis Buffer 2 (5X)*			
	Incubate 30 minutes on an orbital shaker			
Protein Detection	15 $\mu$ l of 4E-BP1 Positive Control Lysate and 15 $\mu$ l of Lysis Buffer 2 (1X) to separate wells			
	Antibody Detection Mix (4X)			
	Cover and incubate 4 hours			
	Read TR-FRET signal			

**Table 5. Assay summary for the 1-plate (all in one well) protocol with adherent or suspension cells**

\*The lysis buffer must be supplemented with the Phosphatase Inhibitor Cocktail (100X) as described in the Reagent Preparation section.

	Untreated Cells	Treated Cells	Positive Control	Negative Control
Suspension cells	8 $\mu$ l	8 $\mu$ l		
Stimulator (3X)		4 $\mu$ l		
Serum-free Medium	4 $\mu$ l			
Positive Control Lysate			15 $\mu$ l	
Lysis Buffer 2 (5X)	3 $\mu$ l	3 $\mu$ l		
Lysis Buffer 2 (1X)				15 $\mu$ l
Antibody Detection Mix (4X)	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l
<b>Total assay volume</b>	20 $\mu$ l	20 $\mu$ l	20 $\mu$ l	20 $\mu$ l

**Table 6. Summary of pipetting protocol 1-Plate (all in one well) protocol using a stimulator**

	Untreated Cells	Treated Cells	Positive Control	Negative Control
Suspension cells	8 $\mu$ l	8 $\mu$ l		
Inhibitor (6X)		2 $\mu$ l		
Stimulator (6X)		2 $\mu$ l		
Serum-free Medium	4 $\mu$ l			
Positive Control Lysate			15 $\mu$ l	
Lysis Buffer 2 (5X)	3 $\mu$ l	3 $\mu$ l		
Lysis Buffer 2 (1X)				15 $\mu$ l
Antibody Detection Mix (4X)	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l	5 $\mu$ l
<b>Total assay volume</b>	20 $\mu$ l	20 $\mu$ l	20 $\mu$ l	20 $\mu$ l

**Table 7. Summary of pipetting protocol 1-Plate (all in one well) protocol using an inhibitor**



## Calculations

1. TR-FRET data are typically calculated and presented ratiometrically using the following formula:

$$[(665 \text{ nm}/615 \text{ nm}) \times 1,000]$$


2. Calculate the TR-FRET ratio for each well.
3. TR-FRET assays are homogeneous; do not subtract average negative control data (no lysate) from any other well readings.
4. For concentration-response curves, analyze data according to a nonlinear regression using the four-parameter logistic equation (sigmoidal dose-response curve with variable slope) and a  $1/Y^2$  data weighting.
5. Assay quality control: The undiluted 4E-BP1 Positive Control Lysate must generate an S/B ratio of at least 2 when compared to the negative control (Lysis Buffer 2 (1X) only). If this is not the case, your reader is not compatible with THUNDER™ TR-FRET Cell Signaling Assay Kits.

*NOTE: The positive control lysate is provided as a control reagent, not for conducting a standard curve.*

## Performance Characteristics

### Representative Data

Data shown here are examples of data typically generated with the 4E-BP1 (Phospho-Thr<sup>37/46</sup>) TR-FRET Assay Kit. The TR-FRET signal was recorded at 665 and 615 nm (EnVision®; lamp excitation) using the recommended plate reader settings. Note that both the TR-FRET ratios and S/B ratios will vary from one TR-FRET-compatible reader to another. In addition, note that excitation with a laser (337 nm) generates higher counts and, usually, higher S/B ratios.

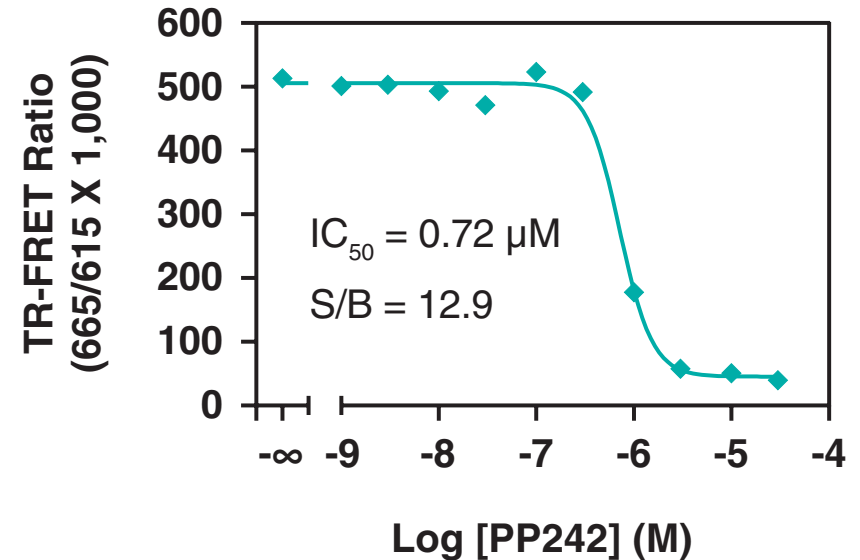

**Figure 4.** Inhibition of 4E-BP1 phosphorylation at Thr<sup>37/46</sup> in A431 cells. A431 cells were seeded at 25,000 cells/well in triplicate and incubated with serial dilutions of the mTOR inhibitor PP242 for 3 hours at 37°C.

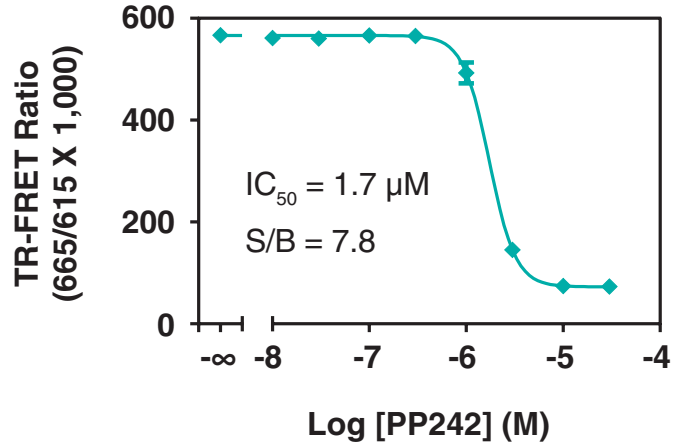


Figure 5. Inhibition of 4E-BP1 phosphorylation Thr<sup>37/46</sup> in HEK293 cells HEK293 cells were seeded at 60,000 cells/well in triplicate and incubated with serial dilutions of PP242 for 3 hours at 37°C.


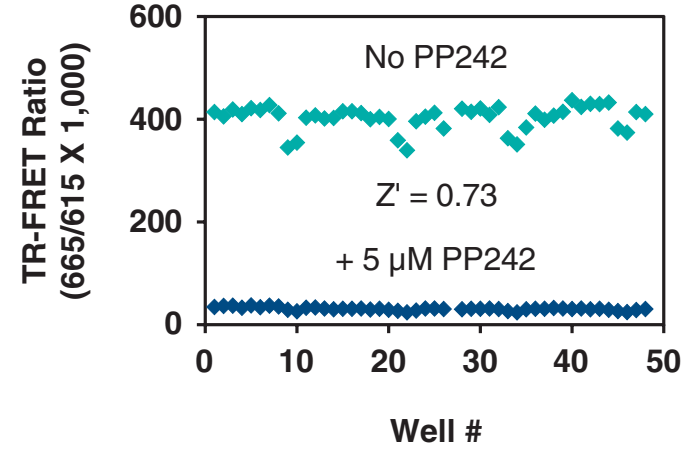


Figure 6. Z'-Factor determination in A431 cells A431 cells were seeded at 25,000 cells/well and incubated without or with 5 μM of PP242 for 3 hours at 37°C. The Z'-factor value was determined using a total of 48 wells for each group. The Z'-factor value of 0.73 indicates that the assay is robust and suitable for high-throughput screening (HTS).

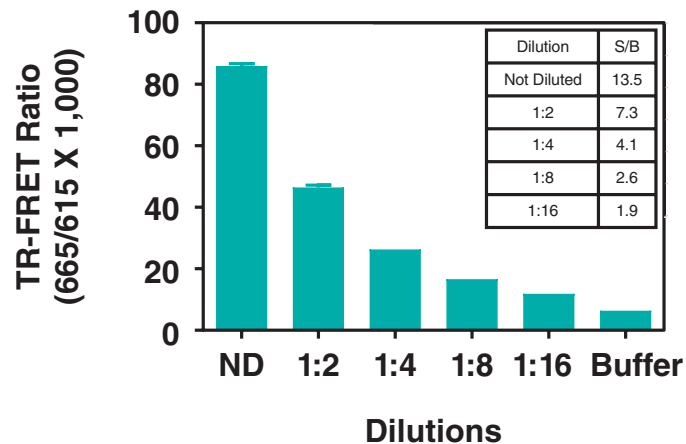

**Figure 7. A431 cell lysate titration (QC Test) for 4E-BP1 (phospho-Thr<sup>37/46</sup>)**  
 The 4E-BP1 (Phospho-Thr<sup>37/46</sup>) TR-FRET Assay Kit is routinely quality control tested using A431 cell lysates. A431 cells were cultured in a T175 flask to 70% confluence. Following cell lysis using 4 ml of Lysis Buffer 2 (1X) lysates were serially diluted with Lysis Buffer 2 (1X) and tested in triplicate.

## RESOURCES

### Troubleshooting

Problem	Possible Causes	Recommended Solutions
Assay S/B ratio is <2 for the positive control lysate versus the negative control (i.e. Lysis Buffer (1X) alone)	Plate reader and/or settings not suitable for TR-FRET assays	Use a filter-based instrument to read the plate(s).
	Use of low-quality water for reagent preparation	Ensure the correct excitation and emission filters and mirror module have been used.
	Use of black plates	Use recommended instrument settings. Optimize the delay time, measurement window, and number of flashes.
	Plate read with the adhesive plate sealer	Only use ultrapure water for preparation of the Lysis and Detection Buffers.  Only use white plates.  The plate sealer MUST be removed before reading the plate(s).

Problem	Possible Causes	Recommended Solutions
Low S/B ratio in the cellular experiment	<p>Suboptimal cell culture and/or treatment conditions</p> <p>Use of a different lysis buffer than the one included in the kit</p> <p>Lack of phosphatase inhibitors in the lysis buffer</p> <p>Use of low-quality water for reagent preparation</p> <p>Use of black plates</p>	<p>Use the positive control lysate to determine whether the poor signal comes from the kit reagents or from the cellular experimental conditions used in the assay.</p> <p>Optimize cell culture conditions. Too high OR low cell numbers can affect basal and maximal activation.</p> <p>Ensure the cell passage number is not too high OR low and that cells are behaving as expected (<i>i.e.</i> doubling time, viability).</p> <p>The lysis buffer MUST be supplemented with the Phosphatase Inhibitor Cocktail (100X) (final concentrations depend on the use of 1X or 5X lysis buffer). Additional phosphatase inhibitors and/or protease inhibitors are typically NOT required.</p> <p>The assay S/B ratio might be increased by decreasing the volume of lysis buffer used to lyse the cells to 25 µl to increase the target protein concentration in the lysate.</p> <p>Only use ultrapure water.</p> <p>Only use white, opaque plates.</p>

## References

1. Musa, J., Orth, M.F., Dallmayer, M., *et al.* Eukaryotic initiation factor 4E-binding protein 1 (4E-BP1): A master regulator of mRNA translation involved in tumorigenesis. *Oncogene* **35(36)**, 4675-4688 (2016).
2. Qin, X., Jiang, B., and Zhang, Y. 4E-BP1, a multifactor regulated multifunctional protein. *Cell Cycle* **15(6)**, 781-786 (2016).
3. Karlsson, E., Pérez-Tenorio, G., Amin, R.P., *et al.* The mTOR effectors 4EBP1 and S6K2 are frequently coexpressed, and associated with a poor prognosis and endocrine resistance in breast cancer: A retrospective study including patients from the randomised Stockholm tamoxifen trials. *Breast Cancer Res.* **15(5)**, R996 (2013).
4. Yeh, C.-J., Chuang, W.-Y., Chao, Y.-K., *et al.* High expression of phosphorylated 4E-binding protein 1 is an adverse prognostic factor in esophageal squamous cell carcinoma. *Virchows Arch.* **458(2)**, 171-178 (2011).
5. Murton, A.J., Alamdari, N., Gardiner, S.M., *et al.* Effects of endotoxaemia on protein metabolism in rat fast-twitch skeletal muscle and myocardium. *PLoS One* **4(9)**, e6945 (2009).
6. Møller, A.B., Lønbro, S., Farup, J., *et al.* Molecular and cellular adaptations to exercise training in skeletal muscle from cancer patients treated with chemotherapy. *J. Cancer Res. Clin. Oncol.* **145(6)**, 1449-1460 (2019).

## NOTES

### Warranty and Limitation of Remedy

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

This document is copyrighted. All rights are reserved. This document may not, in whole or part, be copied, photocopied, reproduced, translated, or reduced to any electronic medium or machine-readable form without prior consent, in writing, from Cayman Chemical Company.

©11/01/2021, Cayman Chemical Company, Ann Arbor, MI, All rights reserved.  
Printed in U.S.A.

