

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com


Product Information



Leukotriene C_4 -d₅ methyl ester

Item No. 9001287

Formal Name:	5S-hydroxy-6R-(S-glutathionyl)-7E,9E,11Z,14Z-d ₅ - eicosatetraenoic acid, methyl ester
Synonym:	LTC ₄ -d ₅ methyl ester
MF:	$C_{31}\dot{H}_{44}D_5N_3O_9S$
FW:	644.8
Chemical Purity:	≥97% LTC ₄ methyl ester
Deuterium	-
Incorporation:	≥99% deuterated forms (d ₁ -d ₅), ≤1% d ₀
Stability:	≥1 year at -80°C
Supplied as:	A solution in ethanol:water (95:5)
UV/Vis.:	λ _{max} : 282 nm
Miscellaneous:	Light Sensitive

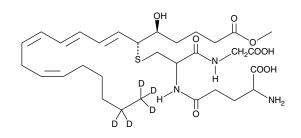

Laboratory Procedures

Leukotriene C_4 - d_5 (LTC₄- d_5) methyl ester contains five deuterium atoms at the 19, 19', 20, 20, and 20 positions. It is intended for use as an internal standard for the quantification of LTC4 methyl ester by GC- or LC-mass spectrometry (MS). For long term storage, we suggest that LTC_4 -d₅ methyl ester be stored as supplied at -80°C. It should be stable for at least one year.

 LTC_4 -d₅ methyl ester is supplied as a solution in ethanol:water (95:5). To change the solvent, simply evaporate the ethanol:water under a gentle stream of nitrogen and immediately add the solvent of choice. Solvents such as DMSO and dimethyl formamide purged with an inert gas can be used. The solubility of LTC_4 -d₅ methyl ester in these solvents is approximately 50 mg/ml.

 LTC_4 -d₅ methyl ester is used as an internal standard for the quantification of LTC_4 methyl ester by stable isotope dilution MS. The accuracy of the sample weight in this vial is between 5% over and 2% under the amount shown on the vial. If better precision is required, the deuterated standard should be quantitated against a more precisely weighed unlabeled standard by constructing a standard curve of peak intensity ratios (deuterated versus unlabeled).

 LTC_4 is the parent cysteinyl-leukotriene produced by the LTC_4 synthase-catalyzed conjugation of glutathione to LTA_4 . LTC₄ is produced by neutrophils, macrophages, and mast cells, and by transcellular metabolism in platelets.¹ It is one of the constituents of slow-reacting substance of anaphylaxis and exhibits potent smooth muscle contracting activity.² LTC₄induced bronchoconstriction and enhanced vascular permeability contribute to the pathogenesis of asthma and acute allergic hypersensitivity.^{3,4} The concentration of LTC₄ required to produce marked contractions of lung parenchymal strips and isolated tracheal rings is about 1 nM.⁴ LTC₄ methyl ester is a more lipid soluble form of LTC₄. The biological activity of LTC₄ methyl ester has not been reported.

References

- 1. Maclouf, J.A. and Murphy, R.C. Transcellular metabolism of neutrophil-derived leukotriene A4 by human platelets. A potential cellular source of leukotriene C₄. J. Biol. Chem. 263, 174-181 (1988).
- 2. Piper, P.J. Formation and actions of leukotrienes. Physiol. Rev. 64, 744-761 (1984).
- 3. Samuelsson, B. Leukotrienes: Mediators of immediate hypersensitivity reactions and inflammation. Science 220, 568-575 (1983).
- 4. Lefer, A.M. Leukotrienes as mediators of ischemia and shock. Biochem. Pharmacol. 35, 123-127 (1986).

Related Products

For a list of related products please visit: www.caymanchem.com/catalog/9001287

WARNING: THIS PRODUCT IS FOR LABORATORY RESEARCH ONLY: NOT FOR ADMINISTRATION TO HUMANS. NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

MATERIAL SAFETY DATA

This material should be considered hazardous until information to the contrary becomes available. Do not ingest, swallow, or inhale. Do not get in eyes, on skin, or on clothing. Wash thoroughly after handling. This information contains some, but not all, of the information required for the safe and proper use of this material. Before use, the user must review the complete Material Safety Data Sheet, which has been sent *via* email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Cayman Chemical Company makes no warranty or guarantee of any kind, whether written or oral, expressed or implied, including without limitation, any warranty of fitness for a particular ose, suitability and merchantability, which extends beyond the description of the chemicals hereof. Cayman warrants only to the original customer that the material will meet our specifications at the time of delivery.

In time of delivery. Cayman will carry out its delivery obligations with due care and skill. Thus, in no event will Cayman have **any obligation or liability**, whether in tort (including negligence) or in contract, for direct, indirect, incidental or consequential damages, even if Cayman is informed about their possible existence. This limitation of liability does not apply in the case of intentional acts or negligence of Cayman, its directors or its employees. Buyer's **exclusive remedy** and Cayman's sole liability hereunder shall be limited to a <u>refund</u> of the purchase price, or at Cayman's option, the <u>replacement</u>, at no cost to Buyer, of all material that

does not meet our specifications. s not need on spectrations. Said refund or preplacement is conditioned on Buyer giving written notice to Cayman within thirty (30) days after arrival of the material at its destination. Failure of Buyer to give said notice within (30) days shall constitute a waiver by Buyer of all claims hereunder with respect to said material. For further details, please refer to our Warranty and Limitation of Remedy located on our website and in our catalog. Copyright Cayman Chemical Company, 11/30/2012

Cayman Chemical

Mailing address

1180 E. Ellsworth Road Ann Arbor, MI 48108 USA

Phone (800) 364-9897 (734) 971-3335

Fax (734) 971-3640

E-Mail

custserv@caymanchem.com

Web

www.caymanchem.com