



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

TECHNICALLY *Speaking*

Place your order with CEDARLANE® or your local distributor.

Please contact CEDARLANE® for lot specific information.

FITC Anti-Rat Endothelium Monoclonal Antibody

CL043F
CL043F-5
LOT: 4331

DESCRIPTION:

Cedarlane's anti-rat endothelium monoclonal antibody recognizes a surface protein of MW 90 kDa and generally reacts with all vascular endothelium in the rat except that of brain capillaries (1). This is the reciprocal tissue pattern to that of the transferrin receptor. It has been shown that the expression of CL043F is on the luminal surface of blood vessels. CL043F labels all peritoneal macrophages, a sub-population of alveolar macrophages (65%) and rare interstitial cells in the brain and heart (1). In addition, Cedarlane's anti-rat endothelium mAb labels circulating erythrocytes, 22% of peripheral blood mononuclear cells and 17% of nucleated cells in bone marrow (1).

CL043F does not label granulocytes, dendritic cells, lymphocytes, or lymphocyte blasts, thymocytes, lymph node cells, mast cells and platelets. CL043F has been invaluable in the demonstration of molecular heterogeneity of vascular endothelium.

PRESENTATION:

100 µg (CL043F) or 500 µg (CL043F-5) FITC conjugated Ig buffered in PBS , 0.02% NaN₃ and EIA grade BSA as a stabilizing protein to bring total protein concentration to 4-5 mg/ml.

STORAGE/STABILITY:

Store at 4°C. For long term storage, aliquot and freeze unused portion at -20°C in volumes appropriate for single usage. Avoid freeze/thaw cycles. Avoid prolonged exposure to light.

For more information or to place an order please contact...

CEDARLANE®
LABORATORIES LIMITED



toll free: 1-800-268-5058
in North America

phone: (905) 878-8891 • fax: (905) 878-7800

5516 - 8th Line, R.R.#2, Hornby, Ontario, CANADA L0P 1E0

or visit our website for a list of our international distributors including contact information

website: www.cedarlanelabs.com • e-mail: info@cedarlanelabs.com

SPECIFICATIONS:

Clone: MRC OX-43

Hybridoma Production:

Immunization: Immunogen: Rat Peritoneal Macrophages
Immunocyte Donor: BALB/c Spleen

Fusion Partner: NSO/U

Specificity: Rat Vascular Endothelium

Isotype: Mouse IgG₁

Presentation: FITC conjugated Ig buffered in PBS, 0.02% NaN₃ and EIA grade BSA as a stabilizing protein to bring total protein concentration to 4-5 mg/ml (Purified from ascitic fluid via Protein G Chromatography).

Antibody Concentration: 0.1 mg/ml

FLOW CYTOMETRY ANALYSIS:**Method:**

1. Prepare a cell suspension in media A. For cell preparations, deplete the red blood cell population with Lympholyte[®]-Rat cell separation medium (CL5040).
2. Wash 2 times.
3. Resuspend the cells to a concentration of 2×10^7 cells/ml in media A. Add 50 μ l of this suspension to each tube (each tube will then contain 1×10^6 cells, representing 1 test).
4. To each tube, add 0.2-0.5 μ g* of **CL043F or CL043F-5**.
5. Vortex the tubes to ensure thorough mixing of antibody and cells.
6. Incubate the tubes for 30 minutes at 4°C.
(It is recommended that the tubes are protected from light, since most fluorochromes are light sensitive.)
7. Wash 2 times at 4°C.
8. Resuspend the cell pellet in 50 μ l ice cold media B.
9. Transfer to suitable tubes for flow cytometric analysis containing 15 μ l of propidium iodide at 0.5 mg/ml in PBS. This stains dead cells by intercalating in DNA.

Media:

- A. Phosphate buffered saline (pH 7.2) + 5% normal serum of host species + sodium azide (100 μ l of 2M sodium azide in 100 mls).
- B. Phosphate buffered saline (pH 7.2) + 0.5% Bovine serum albumin + sodium azide (100 μ l of 2M sodium azide in 100 mls).

Results:Tissue Distribution by Flow Cytometry Analysis:

Rat Strain: Fischer

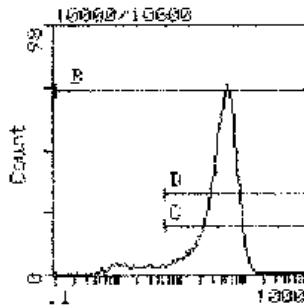
Cell Concentration : 1×10^6 cells per test

Antibody Concentration Used: 0.5 μ g/ 10^6 cells

Isotypic Control: FITC Mouse IgG₁, κ

Cell SourcePercentage of cells stained above control:

| | |
|------------------------|-------|
| Thymus | 3.0% |
| Peritoneal Macrophages | 92.4% |

**LFL1**

Cell Source: Peritoneal Macrophages

Percentage of cells stained above control: 92.4%

N.B. Appropriate control samples should always be included in any labelling studies.

*** For optimal results in various applications, it is recommended that each investigator determine dilutions appropriate for individual use.**

Strain Distribution by Flow Cytometry Analysis:

Procedure: see page 2

Antibody Concentration Used: 0.5 $\mu\text{g}/10^6$ cells

Strains Tested: Wistar, Brown Norway, Buffalo, Fischer 344

Positive: Buffalo, Brown Norway, Fischer 344, Wistar

Negative: none

REFERENCES:

- 1) Robinson, A.P., White, T.M. and D.W. Mason. (1985) Immunology. 57, 231-237. MRC OX-43: a monoclonal antibody which reacts with all vascular endothelium in the rat except that of brain capillaries.

FOR RESEARCH USE ONLY

® is a Registered Trademark of Cedarlane Laboratories Limited.