



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Technically  
Speaking

**CEDARLANE**<sup>®</sup>  
[www.cedarlanelabs.com](http://www.cedarlanelabs.com)



Conveniently Delivering You Today's Innovations  
for the Science of Tomorrow™

**Anti-Human CD44  
Monoclonal Antibody**

Catalogue#	Format	Size	Concentration	Isotype Control
<b>CL7665AP</b>	Purified	250 µg	1.0 mg/ml	CLCMG2A00
<b>CL7665APC</b>	APC	100 µg	0.1 mg/ml	CLCMG2A05
<b>CL7665B</b>	Biotin	100 µg	0.1 mg/ml	CLCMG2A15
<b>CL7665F</b>	FITC	100 µg	0.1 mg/ml	CLCMG2A01
<b>CL7665PE</b>	PE	50 µg	0.1 mg/ml	CLCMG2A04

Isotype: Mouse IgG<sub>2a</sub>

**DESCRIPTION:**

The CD44 antigen is a cell-surface glycoprotein involved in cell to cell and cell-matrix interactions, cell adhesion and migration, which it achieves through its affinity for hyaluronic acid (HA) and possibly also through its affinity for other ligands such as osteopontin, collagens, and matrix metalloproteinases. CD44 is also involved in lymphocyte activation, recirculation and homing, and in hematopoiesis. Altered expression or dysfunction of CD44 is observed in numerous pathogenic phenotypes including various forms of cancer.

**PRESENTATION:**

**Purified:** Purified IgG buffered in PBS and 0.02% NaN<sub>3</sub>. (Purified from ascitic fluid via Protein G Chromatography). For maximum recovery of contents, spin down tube before use.

**APC, Biotin, FITC and PE:** APC/Biotin/FITC/PE conjugated IgG buffered in PBS, 0.02% NaN<sub>3</sub> and EIA grade BSA as a stabilizing protein to bring total protein concentration to 4-5 mg/ml.

**STORAGE/STABILITY:**

For all formats, store at 4°C. DO NOT FREEZE APC and PE conjugates. For long term storage (**Purified, Biotin** and **FITC**), aliquot and freeze unused portion at -20°C in volumes appropriate for single usage. Avoid freeze/thaw cycles.

**APPLICATIONS:**

This antibody is suitable for use in Flow Cytometry (0.2 – 0.1 µg/10<sup>6</sup> cells), Immunohistochemistry, Immunofluorescence, Immunoprecipitation, Western Blotting and Functional Assays (inhibition of CD44-hyaluronate binding).

*Continued.....*

Visit our website for your local distributor.

**CEDARLANE**<sup>®</sup>

[www.cedarlanelabs.com](http://www.cedarlanelabs.com)

An ISO 9001:2015 and ISO 13485:2016  
registered company.

In CANADA: **Toll Free: 1-800-268-5058**

4410 Paletta Court, Burlington, ON L7L 5R2 ph: (289) 288-0001, fax: (289) 288-0020  
e-mail: [general@cedarlanelabs.com](mailto:general@cedarlanelabs.com)

In the USA: **Toll Free: 1-800-721-1644**

1210 Turrentine Street, Burlington, NC 27215 ph: (336) 513-5135, fax: (336) 513-5138  
e-mail: [service@cedarlanelabs.com](mailto:service@cedarlanelabs.com)

**SPECIFICATIONS:**

Clone: BU75

Immunogen: Follicular dendritic cells and germinal center B cells.

Specificity: This antibody is specific for human CD44 (Pgp-1).

**TEST RESULTS:**

Tissue Distribution by Flow Cytometry Analysis:

Cell Concentration:  $1 \times 10^6$  cells per tests


Antibody Concentration Used:  $0.1 \mu\text{g}/10^6$  cells

Cell Source:

Human PBLs

Percentage of cells stained above control:

95.6 %

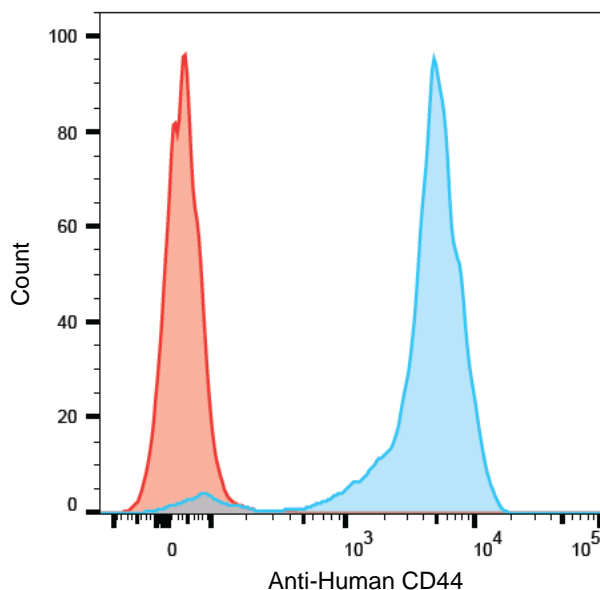

Human peripheral blood lymphocytes were stained with anti-CD44 (clone: BU75) (blue histogram) or Mouse IgG2a isotype control (red histogram).

**N.B. Appropriate control samples should always be included in any labeling studies.**

**\* For optimal results in various applications, it is recommended that each investigator determine dilutions appropriate for individual use.**

**REFERENCES:**

1. Denning, S. M., et al. (1995) in Leucocyte Typing V (Schlossman, S., ed) pp. 1713–1719, Oxford University Press, Oxford.
2. Telen MJ, et al. (1996) A blood group-related polymorphism of CD44 abolishes a hyaluronan-binding consensus sequence without preventing hyaluronan binding. *J Biol Chem.* 271(12):7147-53.
3. Zittermann SI, et al. (2001) Modulation of CD44 in acute lymphoblastic leukemia identifies functional and phenotypic differences of human B cell precursors. *Eur J Haematol.* 66(6):377-82.

**FOR RESEARCH USE ONLY**

® is a registered trademark of Cedarlane