

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in











Mouse Vascular endothelial cell growth factor receptor 2, VEGFR-2/Flk-1 ELISA kit

| Product Code | CSB-E04764m |
|---------------------------------|---|
| Abbreviation | VEGFR-2/Flk-1 |
| Protein Biological Process 1 | Angiogenesis |
| Target Name | kinase insert domain receptor (a type III receptor tyrosine kinase) |
| Uniprot No. | P35918 |
| Alias | CD309, FLK1, VEGFR, VEGFR2, fetal liver kinase-1 kinase insert domain receptor protein-tyrosine kinase receptor Flk-1 soluble VEGFR2 tyrosine kinase growth factor receptor vascular endothelial growt |
| Product Type | ELISA Kit |
| Immunogen Species | Mus musculus (Mouse) |
| Protein Biological Process 3 | Angiogenesis |
| Sample Types | serum, plasma |
| Detection Range | 7.8 pg/mL-500 pg/mL |
| Sensitivity | 1.95 pg/mL |
| Assay Time | 1-5h |
| Sample Volume | 50-100ul |
| Detection Wavelength | 450 nm |
| Lead Time | 3-5 working days |
| Research Area | Cancer |
| Gene Names | Kdr |
| Tag Info | quantitative |
| Protein Description | Sandwich |
| Description | This Mouse VEGFR-2/Flk-1 ELISA Kit was designed for the quantitative measurement of Mouse VEGFR-2/Flk-1 protein in serum, plasma. It is a Sandwich ELISA kit, its detection range is 7.8 pg/mL-500 pg/mL and the sensitivity is 1.95 pg/mL. |
| Target Details | Vascular endothelial growth factor (VEGF) is a major growth factor for endothelial cells. This gene encodes one of the two receptors of the VEGF. This receptor, known as kinase insert domain receptor, is a type III receptor tyrosine |

kinase. It functions as the main mediator of VEGF-induced endothelial

CUSABIO TECHNOLOGY LLC









proliferation, survival, migration, tubular morphogenesis and sprouting. The signalling and trafficking of this receptor are regulated by multiple factors, including Rab GTPase, P2Y purine nucleotide receptor, integrin alphaVbeta3, Tcell protein tyrosine phosphatase, etc.. Mutations of this gene are implicated in infantile capillary hemangiomas.

Product Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

Linearity

To assess the linearity of the assay, samples were spiked with high concentrations of mouse VEGFR-2/Flk-1 in various matrices and diluted with the Sample Diluent to produce samples with values within the dynamic range of the assay.

| ? | Sample | Serum(n=4) |
|---------|-----------|------------|
| 1:100 | Average % | 86 |
| | Range % | 82-91 |
| 1,200 | Average % | 101 |
| 1:200 F | Range % | 98-105 |
| 1:400 | Average % | 87 |
| | Range % | 84-91 |
| 1:800 | Average % | 100 |
| | Range % | 95-105 |

Recovery

The recovery of mouse VEGFR-2/Flk-1 spiked to levels throughout the range of the assay in various matrices was evaluated. Samples were diluted prior to assay as directed in the Sample Preparation section.

| Sample Type | Average % Recovery | Range |
|-------------------|--------------------|---------|
| Serum (n=5) | 89 | 85-94 |
| EDTA plasma (n=4) | 107 | 103-110 |

Typical


These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.









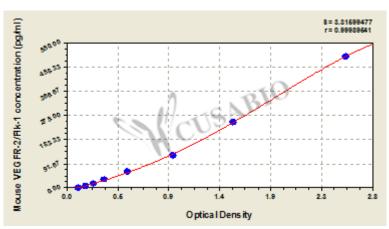

0.063

pg/ml OD1 OD2 Average Corrected

500 2.564 2.536 2.550 2.438 250 1.513 1.531 1.522 1.410 125 0.985 0.961 0.973 0.861 62.5 0.551 0.560 0.556 0.444 31.2 0.345 0.351 0.348 0.236 15.6 0.242 0.251 0.247 0.135

0 0.113 0.110 0.112 ?

7.8 0.174 0.176 0.175