



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



### PRECISION

#### **Intra-assay Precision (Precision within an assay): CV%<15%**

Three samples of known concentration were tested twenty times on one plate to assess.


#### **Inter-assay Precision (Precision between assays): CV%<15%**

Three samples of known concentration were tested in twenty assays to assess.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	20	20	20	20	20	20
Mean(ng/ml)	1.955	1.964	1.945	1.976	1.950	2.003
SD	0.025	0.019	0.028	0.036	0.036	0.033
CV(%)	2.733	2.119	3.052	3.997	4.011	3.685

### TYPICAL DATA

These standard curves are provided for demonstration only. A standard curve should be generated for each set of samples assayed.

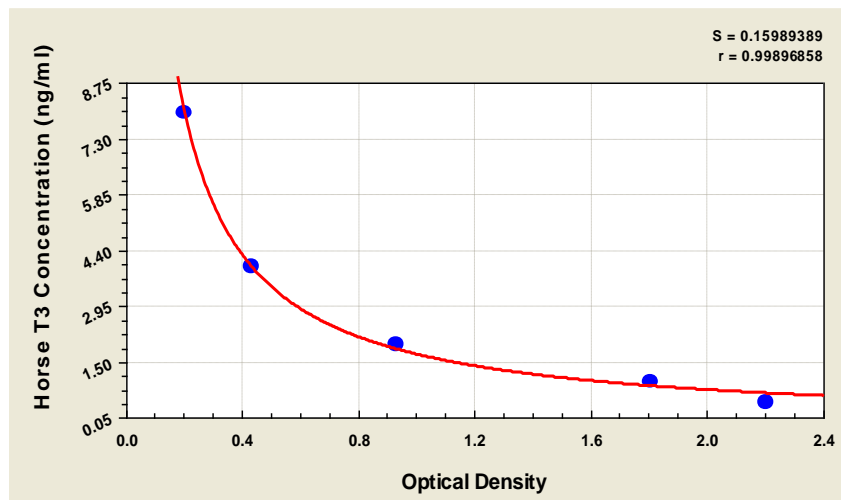

ng/ml	OD1	OD2	Average
0.5	2.198	2.183	2.191
1.0	1.799	1.807	1.803
2.0	0.936	0.958	0.947
4.0	0.467	0.448	0.458
8.0	0.239	0.227	0.233

### LOD

0.25ng/ml