

# Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com




## Sphingosine-1-Phosphate Receptor 3 (S1PR3) ACTOne<sup>™</sup> Stable Cell Line CATALOG NUMBER: CL-01-S1PR3

#### Introduction

S1PR3 is a G-protein-coupled receptor which binds the bioactive signaling molecule sphingosine 1-phosphate (S1P). S1PR3 belongs to a sphingosine-1-phosphate receptor subfamily comprising five members (S1PR3-5). They are expressed in a wide variety of tissues, with each subtype exhibiting a different cell specificity, although they are found at their highest density on leukocytes. S1PR3 has been identified as a functional receptor for sphingosine 1-phosphate and likely contributes to the regulation of angiogenesis and vascular endothelial cell function.

#### Description

Human S1PR3 ACTOne<sup>™</sup> is a HEK-293 CNG cell line that expresses recombinant human S1PR3. HEK-293 CNG cells express a modified CNG (Cyclic Nucleotide Gated) channel that opens in response to elevated intracellular cAMP levels and consequently result in ion flux (often detectable by calcium-responsive dye, <u>Cat# CA-C155</u>) and cell membrane depolarization which can be easily measured with fluorescent Membrane Potential Dye (<u>Cat# CA-M165</u>). The assay allows both end-point and kinetic measurement of intracellular cAMP changes with a FDSS, FLIPR, or a fluorescence microplate reader.

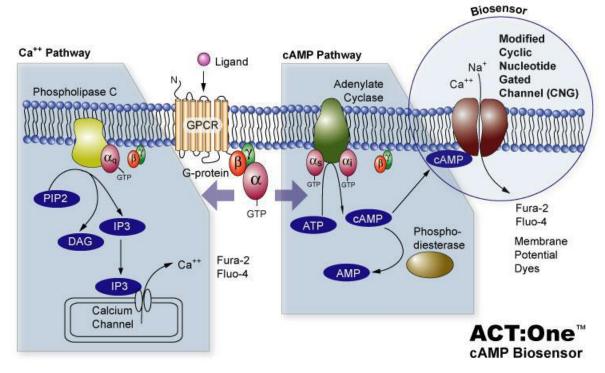

#### **Parental Cells**

HEK-293 CNG cells (originally developed by BD Biosciences by introducing CNG in HEK-293 cells) (Cat# CL-03-PC20)

#### **Gene/Enzyme Introduced**

S1PR3 (Genbank Accession No. NP\_005217.2)

#### Applications

• cAMP dependent human S1PR3 receptor cell based assay



• cell based high-throughput screening of human S1PR3 receptor agonists/antagonists

#### **Functional Test**

- this cell line has been tested positive for S1PR3 receptor specific response
- surviving rate: More than 2.5 million/vial on the second day after thawing
- the receptor specific activity is stable for 10 weeks continuous passage

#### Mycoplasma Contamination Test

This lot of cells has been tested and found to be free of mycoplasma contamination.

#### Content

• Stable S1PR3 receptor cells: 1 mL (1 x 10<sup>6</sup> cells/mL in 70% DMEM, 20% FBS, 10% DMSO)

#### **Growth Properties**

Adherent


#### **Cell Culture Medium**

- Growth medium: DMEM-10% FBS supplemented with 250 µg/ml G418, 1 µg/ml Puromycin
- Freezing medium: 10% DMSO, 90% complete cell culture medium

#### Storage

Remove the frozen cells from the dry ice packaging and immediately place the cells at a temperature below -130°C, preferably in liquid nitrogen vapor, until ready for use.

#### Data Example

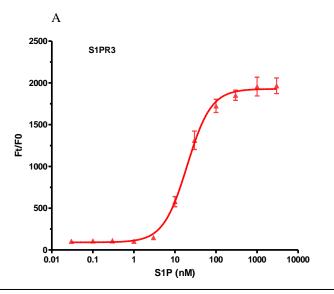

#### Figure 1. Response of ACTOne<sup>™</sup> S1PR3 cell line

ACTOne<sup>™</sup> S1PR3 receptor cells and parental cells (Cat# CL-03-PC20) were plated overnight in 20 µl culture medium on a 384 well Biocoat plate. The next day, cells were dye-loaded with 20 µl/well of 1x Dye-loading solution (Elite<sup>™</sup> Non-Wash Calcium Dye, Cat# CA-C155). After 1 hour of incubation at 37 °C, the plate was read on Hmamastu FDSS 7000 with online addition of the ligand. Fluorescent intensity was used to plot the figure.

A. Dose response curve of S1P in ACTOne<sup>™</sup> S1PR3 cell line. EC50 = 20 nM \*Parental cells do not respond to S1P (data not shown)



#### **Notice to Purchaser**

This cell line is sold with a use license. It may not be transferred to third parties, resold, modified for resale without written approval of eEnzyme LLC. Refer to the license agreements for details on the usage restrictions.