

# Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

# Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

# SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com

# **Elabscience**®

# FITC Anti-Human CD86 Antibody[IT2.2]

Catalog No.E-AB-F1269CStorageStore at 2~8°C, Avoid freeze / thaw cycles

ReactivityHumanApplicationsFCM

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.

### **Antigen Information**

Alternate Names Uniprot ID Gene ID Background	<ul> <li>T-lymphocyte activation antigen CD86,CD28LG2,B7-2,B70</li> <li>P42081</li> <li>942</li> <li>CD86 is an 80 kD immunoglobulin superfamily member also known as B7-2, B70, and Ly-58.</li> <li>CD86 is expressed on activated B and T cells, monocytes/macrophages, dendritic cells, and astrocytes. CD86, along with CD80, is the ligand of CD28 and CD152 (CTLA-4). CD86 is expressed earlier in the immune response than CD80. CD86 has also been shown to be involved in immunoglobulin class-switching and triggering of NK cell-mediated cytotoxicity. CD86 binds</li> </ul>

#### **Product Details**

Form	Liquid
Size	20Tests/100Tests/100Tests×2
Clone No.	IT2.2
Host	Mouse
Isotype	Mouse IgG2b, ĸ
Reactivity	Human
Application	FCM
Isotype Control	FITC Mouse IgG2b, κ Isotype Control[MPC-11] [Product E-AB-F09812C]
Storage Buffer	Phosphate buffered solution, pH 7.2, containing 0.09% stabilizer and 1% protein protectant.
Shipping	Biological ice pack at 4 °C
Stability & Storage	Keep as concentrated solution.
	Store at 2~8°C and protected from prolonged exposure to light.Do not freeze.
	This product is guaranteed up to one year from purchase.

**For Research Use Only** 

# **Elabscience**®

## Fluorophore


#### **Conjugation:** FITC

FITC is designed to be excited by the Blue laser (488 nm) and detected using an optical filter centered near 530 nm (e.g., a 525/40 nm bandpass filter).

### **Recommended usage**

Each lot of this antibody is quality control tested by flow cytometric analysis. The amount of the reagent is suggested to be used 5  $\mu$ L of antibody per test (million cells in 100  $\mu$ L staining volume or per 100  $\mu$ L of whole blood). Please check your vial before the experiment. Since applications vary, the appropriate dilutions must be determined for individual use.

### **Product data**

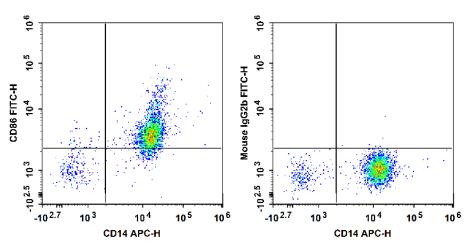

Human peripheral blood are stained with APC Anti-Human CD14 Antibody and FITC Anti-Human CD86 Antibody (Left). Cells in the monocyte gate were used for analysis. Cells are stained with APC Anti-Human CD14 Antibody and FITC Mouse IgG2b, κ Isotype Control (Right).

### **Related Information**

- 1. Sample Preparation for Flow Cytometry https://www.elabscience.com/List-detail-5594.html
- 2. Staining Cell Surface Targets for Flow Cytometry https://www.elabscience.com/List-detail-5568.html
- 3. Flow Cytometry Troubleshooting Tips <u>https://www.elabscience.com/List-detail-5593.html</u>
- 4. How to select the appropriate detection channel through the spectrogram? <u>https://www.elabscience.com/List-detail-459742.html</u>