



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Recombinant Swine Flt-3 Ligand protein(His Tag)



Catalog Number:PKSS000013

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Synonyms	FLT3LG;Fms-related tyrosine kinase 3 ligand transcript variant 1
Species	Porcine
Expression Host	E.coli
Sequence	Ser 30-Leu 205
Accession	D2K7D6
Calculated Molecular Weight	20.4 kDa
Observed molecular weight	17 kDa
Tag	C-His
Bioactivity	Measure by its ability to induce proliferation in OCI-AML5 cells.The ED ₅₀ for this effect is < 5 ng/mL.

Properties

Purity	> 95 % as determined by reducing SDS-PAGE.
Endotoxin	< 0.1 EU per µg of the protein as determined by the LAL method.
Storage	Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of reconstituted samples are stable at < -20°C for 3 months.
Shipping	This product is provided as lyophilized powder which is shipped with ice packs.
Formulation	Lyophilized from sterile PBS, pH 7.4. Normally 5 % - 8 % trehalose, mannitol and 0.01% Tween80 are added as protectants before lyophilization. Please refer to the specific buffer information in the printed manual.
Reconstitution	Please refer to the printed manual for detailed information.

Data

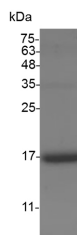

> 95 % as determined by reducing SDS-PAGE.

Background

The cluster of differentiation (CD) system is commonly used as cell markers in immunophenotyping. Different kinds of cells in the immune system can be identified through the surface CD molecules which associating with the immune function of the cell. There are more than 320 CD unique clusters and subclusters have been identified. Some of the CD molecules serve as receptors or ligands important to the cell through initiating a signal cascade which then alter the behavior of the cell. Some CD proteins do not take part in cell signal process but have other functions such as cell adhesion. CD135, also known as FLT-3, FLK-2, is a member of the CD system. CD135 is an important cell surface

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Toll-free: 1-888-852-8623

Tel: 1-832-243-6086

Fax: 1-832-243-6017

Web: www.elabscience.com

Email: techsupport@elabscience.com

Recombinant Swine Flt-3 Ligand protein(His Tag)

Catalog Number:PKSS000013



marker recognized by specific sets of antibodies to identify the types of hematopoietic (blood) progenitors in the bone marrow and its function to differentiate hematopoietic stem cells, which are CD135 negative, from multipotent progenitors, which are CD135 positive. CD135 is a receptor tyrosine kinase type III for the cytokine Flt3 ligand and activates signaling through second messengers by binding to Flt3. Signaling through CD135 is important for lymphocyte development. The encoding gene CD135 is a proto-oncogene to which mutations have occurred will lead to cancer such as leukemia.

For Research Use Only

A Reliable Research Partner in Life Science and Medicine

Toll-free: 1-888-852-8623

Tel: 1-832-243-6086

Fax: 1-832-243-6017

Web: www.elabscience.com

Email: techsupport@elabscience.com