



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



## UK Office

### Everest Biotech Ltd

Cherwell Innovation Centre  
77 Heyford Park  
Upper Heyford  
Oxfordshire  
OX25 5HD  
UK

Enquiries:

[info@everestbiotech.com](mailto:info@everestbiotech.com)

Sales:

[sales@everestbiotech.com](mailto:sales@everestbiotech.com)

Tech support:

[support@everestbiotech.com](mailto:support@everestbiotech.com)

Tel: +44 (0)1869 238326

Fax: +44 (0)1869 238327

## US Office

### Everest Biotech c/o Abcore

405 Maple Street, Suite A106  
Ramona,  
CA 92065  
USA

Inquiries:

[info@everestbiotech.com](mailto:info@everestbiotech.com)

Sales:

[usasales@everestbiotech.com](mailto:usasales@everestbiotech.com)

Tech support:

[support@everestbiotech.com](mailto:support@everestbiotech.com)

Tel: 888-320-4628 (toll-free)

Fax: 888-841-9041

[www.everestbiotech.com](http://www.everestbiotech.com)

**Research Use Only. Not for  
diagnostic or therapeutic use.**

## EB05035 - Goat Anti-Aurora kinase A / AURKA Antibody

Size: 100µg specific antibody in 200µl



### Target Protein

**Principal Names:** AURKA, aurora kinase A, AIK, ARK1, AURA, AURORA2, BTAK, MGC34538, STK15, STK6, STK7, IPL1-related kinase, aurora-A, aurora-related kinase 1, aurora/IPL1-like kinase, breast-tumor-amplified kinase, serine/threonine kinase 15, serine/threonine kinase 6, serine/threonine protein kinase 15, serine/threonine protein kinase 6, OTTHUMP00000031340, OTTHUMP00000031342, OTTHUMP00000031343, OTTHUMP00000031344, OTTHUMP00000031345, OTTHUMP00000166071

**Official Symbol:** AURKA

**Accession Number(s):** NP\_003591.2

**Human GeneID(s):** [6790](#)

**Important Comments:** Reported variants represent identical protein: NP\_003591.2; NP\_940835.1; NP\_940836.1; NP\_940837.1; NP\_940838.1; NP\_940839.1

### Immunogen

Peptide with sequence C-QNKESASKQS, from the C Terminus of the protein sequence according to NP\_003591.2.

Please note the [peptide](#) is available for sale.

### Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

### Applications Tested

**Peptide ELISA:** antibody detection limit dilution 1:128000.

**Western blot:** Approx 45kDa band observed in lysates of cell line HeLa and Jurkat lysate (calculated MW of 45.8kDa according to NP\_003591.2). Recommended concentration: 0.1-0.3µg/ml.

**Immunoprecipitation:** This antibody was deemed fit for IP under native conditions (observations from a customer).

**Additional validation:** This antibody has been successfully used in the following paper: Sikorski et al. (2018) PMID: 30377371.

### Species Reactivity

**Tested:** Human

**Expected from sequence similarity:** Human

### Specific Reference


**This antibody has been successfully used in the following paper:**

Krzysztof Sikorski, Adi Mehta, Marit Inngjerdigen, Flourina Thakor, Simon Kling, Tomas Kalina, Tuula A. Nyman, Maria Ekman Stensland, Wei Zhou, Gustavo A. De Souza, Lars Holden, Jan Stuchly, Markus Templin and Fridtjof Lund-Johansen

A high-throughput pipeline for validation of antibodies

Nat Methods. 2018 Nov;15(11):909-912

PMID: 30377371

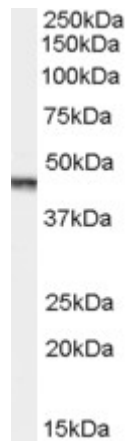

EB05035 (0.1µg/ml) staining of Jurkat lysate (35µg protein in RIPA buffer). Primary incubation was 1 hour.  
Detected by chemiluminescence.