



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Vibrio parahaemolyticus TRH Toxin antibody [vp-01]

Cat No. GTX00855

Host	Mouse
Clonality	Monoclonal
Isotype	IgG1
Application	WB, ELISA
Reactivity	Vibrio parahaemolyticus

Package
100 µg

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Dilution
WB	1:500-1:1000
ELISA	Assay dependent

Not tested in other applications.

Specificity/Sensitivity This antibody is useful for detecting food-poisoning V. hemolyticus strains.


PROPERTIES

Form	Liquid
Buffer	Filter-sterilized PBS, 50% glycerol
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.5 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Culture supernatant of V. parahemolyticus
Purification	Protein A/G purified
Conjugation	Unconjugated
Note	For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.



For full product information, images and publications, please visit our [website](#).

DATA IMAGES

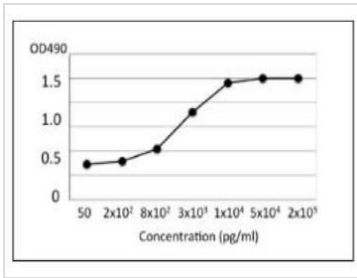


GTX00855 ELISA Image

ELISA analysis of culture medium of *V. parahaemophilus* trh+ using GTX00855 *Vibrio parahaemolyticus* TRH Toxin antibody [vp-01].

Protein coating amount : 100µl per well at indicated amount

Antibody amount : 100µl (1:500) per well

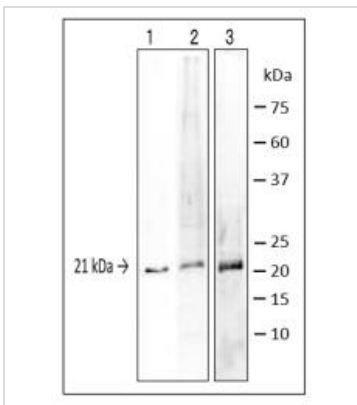


GTX00855 WB Image

WB analysis of culture medium of *V. parahaemolyticus* (trh+) using GTX00855 *Vibrio parahaemolyticus* TRH Toxin antibody [vp-01].

Lane 1-2 : GTX00855

Lane 3 : GTX64192

Dilution : 1:1000

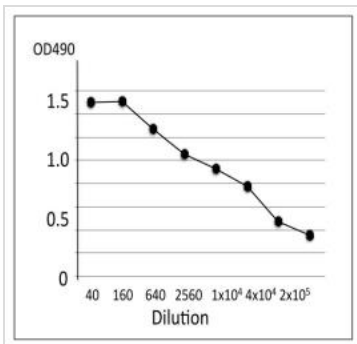

GTX00855 ELISA Image

ELISA analysis of culture medium of *V. parahaemophilus* trh+ using GTX00855 *Vibrio parahaemolyticus* TRH Toxin antibody [vp-01].

Protein coating amount : 100µl (1 µg/ml) per well

Antibody amount : 100µl at indicated dilution

Tabel. Reactivity of clone[vp-01] with various food poisoning bacteria.

	ELISA	WB
<i>Vibrio parahaemolyticus</i> (NBR32711)	+	21K
Other 3 isolated strains	+	21K
<i>Salmonella</i> Enteritidis	-	-
<i>E coli</i> O157:H7	-	-
<i>Staphylococcus aureus</i>	-	-
<i>Bacillus cereus</i>	-	-
Partially purified TDH	+	-

TDH: Thermostable direct hemolysin

GTX00855 Image

Reactivity of GTX00855 *Vibrio parahaemolyticus* TRH Toxin antibody [vp-01] with various food poisoning bacteria.



For full product information, images and publications, please visit our [website](#).