



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Plasmodium falciparum Ferredoxin antibody

Cat No. GTX00920

Host	Rabbit
Clonality	Polyclonal
Isotype	IgG
Application	WB, ICC/IF, ELISA
Reactivity	Plasmodium falciparum

Package
200 µg

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Dilution
WB	1:500-1:2000
ICC/IF	Assay dependent
ELISA	Assay dependent

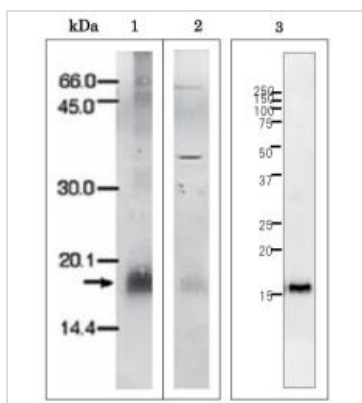
Not tested in other applications.

PROPERTIES

Form	Liquid
Buffer	Filter-sterilized PBS, 50% glycerol
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	4 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Ferredoxin (Pd) protein purified from Malaria parasite, Plasmodium falciparum
Purification	Protein A purified
Conjugation	Unconjugated
Note	For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.



For full product information, images and publications, please visit our [website](#).

DATA IMAGES

GTX00920 WB Image

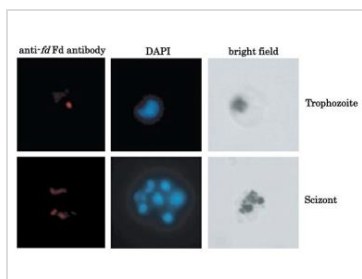
WB analysis of various samples using GTX00920 Plasmodium falciparum Ferredoxin antibody. Molecular mass of Pf Fd is 18 kDa.

Lane 1 : Purified recombinant P. falciparum Ferredoxin (10 ng)

Lane 2 : Partially purified Ferredoxin protein from culture of P. falciparum

Lane 3 : Purified recombinant P. falciparum Ferredoxin (1.4 ng)

Dilution : 1:1000


GTX00920 ICC/IF Image

WB analysis of Trophozoite and schizont stages of P. falciparum using GTX00920 Plasmodium falciparum Ferredoxin antibody. Dark spots in bright field microscopy (left panels) are hemozoin pigment.



For full product information, images and publications, please visit our [website](http://www.genetex.com).