



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Melan A antibody [MSVA-901M+] HistoMAX™

Cat. No. GTX04384

Host	Mouse
Clonality	Monoclonal
Isotype	IgG1
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
--------------------	----------------------

IHC-P	1:50-1:100
-------	------------

Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 9.0 Tris-EDTA-based Target Retrieval Solution buffer. Based on testing result, pH 9.0 showed a strongest staining pattern. However, pH 9.0 results in only a slight reduction of the staining intensity as compared to pH 9.0. Therefore, the pH 7.8 was an optimal for manual staining because of the better tissue preservation at pH 7.8 than at pH 9.0.

This antibody has been validated by Autostainer including Agilent/Dako - Autostainer Link 48, Bond RX-Lecia, and Roche-Ventana Discovery ULTRA. The detailed steps by Autostainer please refer to the protocol.

Not tested in other applications.

Calculated MW	13 kDa. (Note)
---------------	----------------------------------

PROPERTIES

Form	Liquid
------	--------

Buffer	PBS, 0.05% BSA (Please contact us for PBS only format)
--------	--

Preservative	0.05% sodium azide
--------------	--------------------

Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
---------	--

Concentration	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
---------------	--

Immunogen	Recombinant full-length hMART-1 protein
-----------	---



For full product information, images and publications, please visit our [website](#).


Purification Protein A/G purified

Conjugation Unconjugated

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.


DATA IMAGES

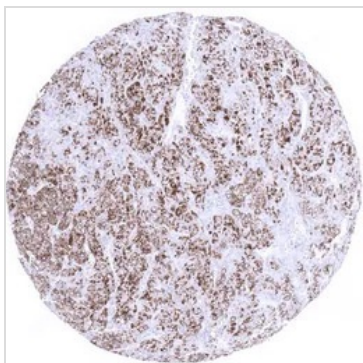

GTX04384 IHC-P Image

IHC-P analysis of human adrenocortical adenoma tissue section using GTX04384 Melan A antibody [MSVA-901M+] HistoMAX™.


Adrenocortical adenoma showing diffuse strong Melan A positivity.

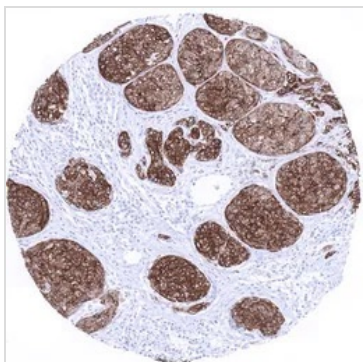

GTX04384 IHC-P Image

IHC-P analysis of human skin Malignant melanoma tissue section using GTX04384 Melan A antibody [MSVA-901M+] HistoMAX™.

Skin Malignant melanoma with strong Melan A positivity in all tumor cells.

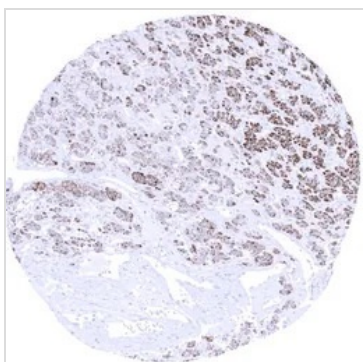

GTX04384 IHC-P Image

IHC-P analysis of human clear cell renal cell carcinoma tissue section using GTX04384 Melan A antibody [MSVA-901M+] HistoMAX™.

Clear cell renal cell carcinoma showing moderate to strong Melan A immunostaining of tumor cells.



For full product information, images and publications, please visit our [website](#).