



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Progesterone Receptor antibody [MSVA-570R] HistoMAX™

Cat. No. GTX04431

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:50-1:100

Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

This antibody has been validated by Autostainer including Agilent/Dako - Autostainer Link 48, Bond RX-Lecia, and Roche-Ventana Discovery ULTRA. The detailed steps by Autostainer please refer to the protocol.

Not tested in other applications.

Calculated MW 99 kDa. ([Note](#))

PROPERTIES

Form	Liquid
Buffer	PBS, 0.05% BSA (Please contact us for PBS only format)
Preservative	0.05% sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant fragment (around aa 1-200) of human PGR
Purification	Protein A/G purified




For full product information, images and publications, please visit our [website](#).

Date 2024 / 06 / 27 Page 1 of 2

Conjugation Unconjugated


Note For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.
Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

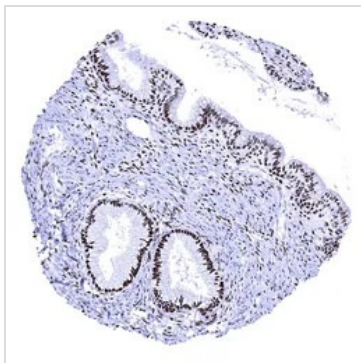


GTX04431 IHC-P Image

IHC-P analysis of human uterus endocervix tissue section using GTX04431 Progesterone Receptor antibody [MSVA-570R] HistoMAX™.
Strong PR immunostaining in epithelial and stromal cells of the endocervix.

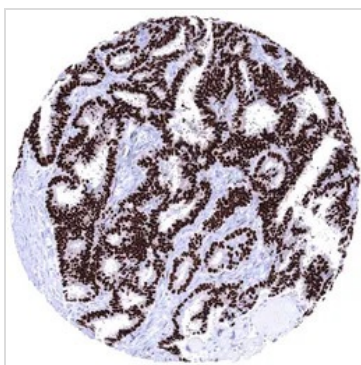

GTX04431 IHC-P Image

IHC-P analysis of human endometrioid ovarian carcinoma tissue section using GTX04431 Progesterone Receptor antibody [MSVA-570R] HistoMAX™.
Endometrioid carcinoma of the ovary depicting a strong PR immunostaining of all tumor cells.

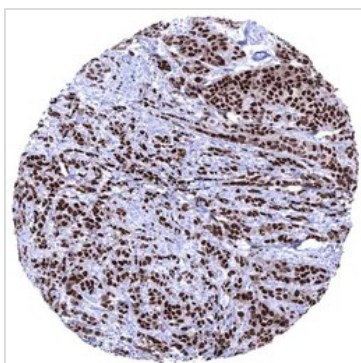

GTX04431 IHC-P Image

IHC-P analysis of human invasive breast cancer tissue section using GTX04431 Progesterone Receptor antibody [MSVA-570R] HistoMAX™.
Invasive breast cancer of no special type (NST) with strong PR immunostaining in 100 of tumor cells.



For full product information, images and publications, please visit our [website](#).