



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

LI-Cadherin antibody [MSVA-517M] HistoMAX™

Cat. No. GTX04447

Host	Mouse
Clonality	Monoclonal
Isotype	IgG2b
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)
[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Manual staining : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Tris-EDTA-based Target Retrieval Solution buffer.

This antibody has been validated by Autostainer including Agilent/Dako - Autostainer Link 48, Bond RX-Lecia, and Roche-Ventana Discovery ULTRA. The detailed steps by Autostainer please refer to the protocol.

Not tested in other applications.

Calculated MW 92 kDa. ([Note](#))

PROPERTIES

Form	Liquid
Buffer	PBS, 0.05% BSA (Please contact us for PBS only format)
Preservative	0.05% sodium azide
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.2 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Recombinant fragment (around aa 242-418) of human Cadherin 17 protein (CDH17)
Purification	Protein A/G purified



For full product information, images and publications, please visit our [website](#).


Date 2024 / 06 / 27 Page 1 of 2

Conjugation Unconjugated

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.


Note Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

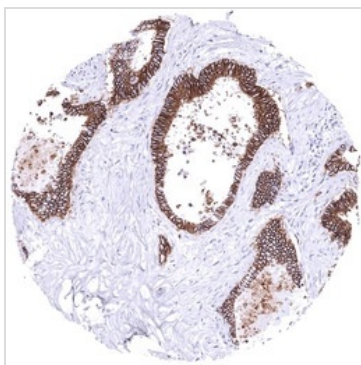

GTX04447 IHC-P Image

IHC-P analysis of human colon colorectal adenocarcinoma tissue section using GTX04447 LI-Cadherin antibody [MSVA-517M] HistoMAX™. Colon Colorectal adenocarcinoma with strong LI-Cadherin immunostaining of all tumor cells LI-Cadherin immunohistochemistry.

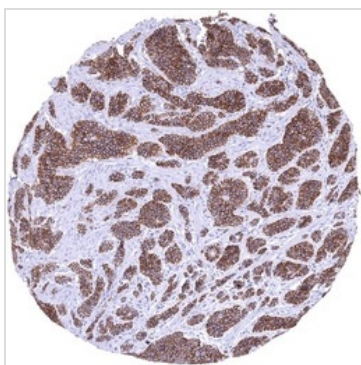

GTX04447 IHC-P Image

IHC-P analysis of human ileum Neuroendocrine tissue section using GTX04447 LI-Cadherin antibody [MSVA-517M] HistoMAX™. Ileum Neuroendocrine tumor with strong LI-Cadherin immunostaining of all tumor cells LI-Cadherin immunohistochemistry.



GTX04447 IHC-P Image

IHC-P analysis of human rectum mucosa tissue section using GTX04447 LI-Cadherin antibody [MSVA-517M] HistoMAX™. A strong membranous LI-Cadherin staining occurs in all epithelial cells of the rectum mucosa LI-Cadherin immunohistochemistry.



For full product information, images and publications, please visit our [website](#).