



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

SignalPlus Antibody Enhancer

Cat No. GTX49999

Application WB, Dot, ELISA

Reference (3)

Package

50 ml

PRODUCT

Product Description

SignalPlus Antibody Enhancer for Western Blotting and ELISA enhances antigen-antibody reactions. SignalPlus Antibody Enhancer works well with chemiluminescence, colorimetric and fluorescence detection system. It can significantly enhance detection of weak immunoreactive and low abundance proteins in a variety of immunoassays such as Western blotting, dot blotting and ELISA. Simply dilute antibodies with SignalPlus Antibody Enhancer and process the rest of the procedures as usual. No additional steps are required. Signal enhancement is protein dependent and could vary from several folds to more than ten-fold.

APPLICATION

Application Note

Dilute primary antibody in Solution 1 and secondary antibody in Solution 2. We recommend using a 10-fold lower working concentration of primary antibody than is usually required for western blot. Optimal working dilutions should be determined experimentally by the end user.


PROPERTIES

Form Liquid

Storage Store at 4°C.

Note For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

DATA IMAGES

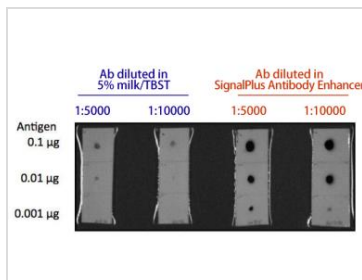


GTX49999 Dot Image

Antibody diluted in SignalPlus Antibody Enhancer solution can detect more than 100-fold lower levels of antigen than antibody diluted in standard buffers




For full product information, images and publications, please visit our [website](#).

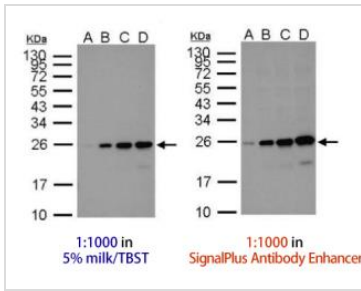

GTX49999 WB Image

Antibody diluted in SignalPlus Antibody Enhancer solution has superior signal to antibody diluted in standard buffers, even at a 10-fold lower concentration of antibody.

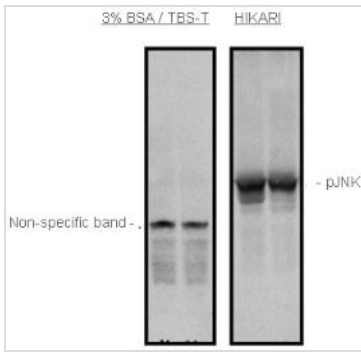

GTX49999 WB Image

Detection enhancement of pJNK with treatment of Signal + in Western blot. Treatment with GTX49999 resulted in clear signal enhancement and background suppression compared with the conventional method using TBS-T. Primary Ab: Anti-pJNK (1: 1,000) Secondary Ab: Anti-rabbit IgG - HRP (1: 5,000) Sample: Mouse Embryonic Fibroblast (MEF) Detection: SuperSignal West Pico



GTX49999 WB Image

Detection enhancement of Rac1 with treatment of Signal + in Western blot. Treatment with GTX49999 resulted in clear signal, whereas dilution with conventional PBS-T resulted in no signal. Primary Ab: Anti-Rac1 (1: 1,000) Secondary Ab: Anti-Mouse IgG-HRP (1: 2,500) Sample: Rat cortical primary neuron Detection: ECL



For full product information, images and publications, please visit our [website](http://www.genetex.com).