



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



Anterior Gradient 2 antibody [HMV325] HistoMAX™

Cat. No. GTX640557

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)

[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Target Retrieval Solution buffer.

Not tested in other applications.

PROPERTIES

Form	Liquid
Buffer	PBS, 2% BSA, 20% Glycerol
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.02 mg/ml
Immunogen	Carrier-protein conjugated synthetic peptide encompassing a sequence within the N-terminus region of human Anterior Gradient 2. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated

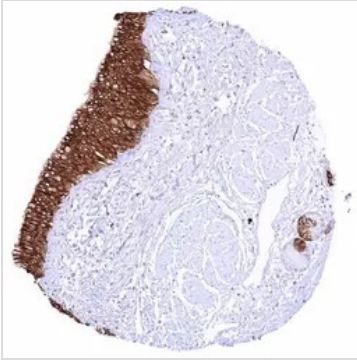


For full product information, images and publications, please visit our [website](#).

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

GTx640557 IHC-P Image

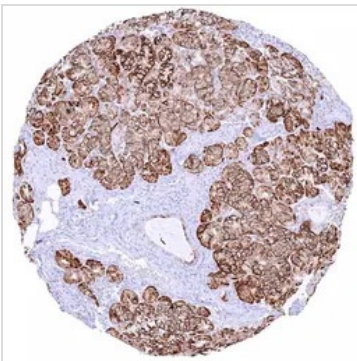
IHC-P analysis of human bronchial mucosa tissue using GTx640557 Anterior Gradient 2 antibody [HmV325] HistoMAX™.

Strong predominantly cytoplasmic anterior gradient 2 staining of all epithelial cells.


GTx640557 IHC-P Image

IHC-P analysis of human stomach tissue using GTx640557 Anterior Gradient 2 antibody [HmV325] HistoMAX™.

Strong predominantly cytoplasmic anterior gradient 2 staining of all epithelial cells except parietal cells which only show a faint staining.


GTx640557 IHC-P Image

IHC-P analysis of human sublingual gland from oral cavity tissue using GTx640557 Anterior Gradient 2 antibody [HmV325] HistoMAX™.

Strong predominantly cytoplasmic anterior gradient 2 staining of subset of mucinous glandular cells. Weak to moderate AGR2 staining of a subset of excretory duct cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).