



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

CA19-9 antibody [HMV333] HistoMAX™

Cat. No. GTX640565

Host	Mouse
Clonality	Monoclonal
Isotype	IgG1
Application	IHC-P
Reactivity	Species independent

Package

500 µl, 100 µl

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Target Retrieval Solution buffer.

Not tested in other applications.

PROPERTIES

Form	Liquid
Buffer	PBS, 2% BSA, 20% Glycerol
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.22 mg/ml
Immunogen	The immunogen used to generate this antibody corresponds to CA19-9.
Purification	Affinity purified by Protein G.
Conjugation	Unconjugated

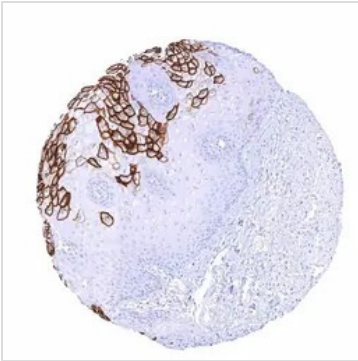
Note

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

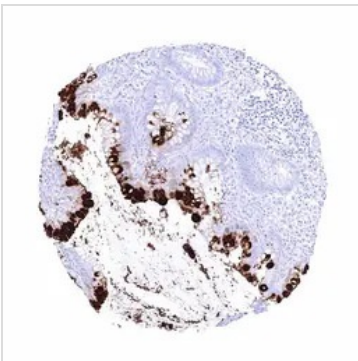
Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.



For full product information, images and publications, please visit our [website](#).

DATA IMAGES

GTX640565 IHC-P Image

IHC-P analysis of human esophagus tissue using GTX640565 CA19-9 antibody [HMOV333] HistoMAX™. Strong membranous CA19-9 staining of a significant fraction of cells in superficial layers of the squamous epithelium.


GTX640565 IHC-P Image

IHC-P analysis of human mucosa from appendix tissue using GTX640565 CA19-9 antibody [HMOV333] HistoMAX™. Moderate to strong apical membranous CA19-9 staining of most superficial epithelial cells. In addition, strong membranous and cytoplasmic positivity of a large fraction of goblet cells.


GTX640565 IHC-P Image

IHC-P analysis of human pancreatic ductal adenocarcinoma (PDAC) tissue using GTX640565 CA19-9 antibody [HMOV333] HistoMAX™. Ductal adenocarcinoma with strong CA19-9 staining of tumor cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).