



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien


T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

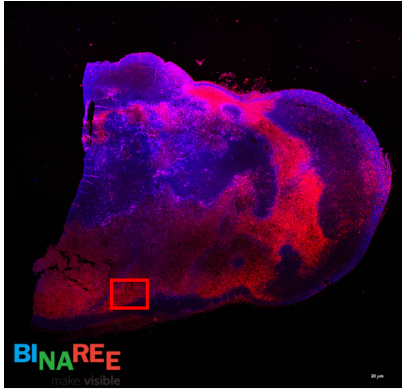

Anti-ICAM-1 In Vivo Antibody – Low Endotoxin (15.2) [ICH1027]

Description

Bulk anti-ICAM-1 antibody clone 15.2

Product Benefits:

ichorbio's anti-ICAM-1 In Vivo Antibody – Low Endotoxin (15.2) [ICH1027] is manufactured in a cGMP compliant, ISO Quality Standard 9001:2015 facility. ichorbio's low endotoxin antibodies have half the endotoxin of comparable antibodies from our [competitors](#) at less than 1.0 EU/mg. If ichorbio's low endotoxin antibodies are not low enough we also offer ultra low endotoxin antibodies which have even less endotoxin (<0.75EU/mg) at an even higher purity (98% versus 95%). ichorbio: the best antibodies for *in vivo* research.

Target:

ICAM-1

Clone:

15.2

Isotype:

Mouse IgG1

Other Names:

Intercellular adhesion molecule 1, ICAM1, Major group rhinovirus receptor, CD54

Host:

Mouse

Species Reactivity:

Human, Mouse, Rat

Specificity:

Mouse anti-ICAM-1 in vivo antibody – low endotoxin (15.2) recognizes an epitope on Human ICAM-1

Purification Method:

This bulk anti-ICAM-1 antibody clone 15.2 was purified using multi-step affinity chromatography methods such as Protein A or G depending on the species and isotype.

Concentration:

0.5 mg/ml

Formulation:

0.01 M phosphate buffered saline (PBS) pH 7.2, 150 mM NaCl with no carrier protein, potassium or preservatives added. BSA and Azide free.

Purity:

>95% by SDS-PAGE and HPLC

>98% by SDS-PAGE and HPLC

Endotoxin:

? 1.0 EU/mg as determined by the LAL method

? 0.75 EU/mg as determined by the LAL method

Aggregation:

Aggregation level ? 5%

Aggregation level ? 1%

Storage:

Anti-ICAM-1 In Vivo Antibody – Low Endotoxin (15.2) is stable for at least one week when stored at 2-8°C. For long term storage, aliquot in working volumes without diluting and store at – 20°C in a manual defrost freezer. Avoid Repeated Freeze Thaw Cycles.

Applications:

IHC (FFPE), Flow Cytometry, Immunoprecipitation, Western Blot

Application Notes:

Each investigator should determine their own optimal working dilution for specific applications.

Use:

Products are for research use only. Not for use in diagnostic or therapeutic procedures.

Antibodies against the same target:

[Anti-ICAM-1 In Vivo Antibody – Low Endotoxin \[BE29G1\] \(ICH1101\)](#)

Immunofluorescence (paraffin-embedded sections):

Immunofluorescence analysis of paraffin-embedded breast cancer tissue section labeling ICAM-1 (1:100) overnight at 4°C, followed by goat anti-mouse IgG H&L (Alexa Fluor® 647-red) secondary antibody (1:500 dilution). MDA-MB-231 cells were injected via mouse nipple and grown for 5 weeks. The nuclear counter stain is DAPI (blue). Image was acquired on a Nikon A1R microscope system at 4x magnification (first image) or 60x magnification (second image).