



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

4-Nitrophenyl phosphate disodium hexahydrate

Cat. No.:	HY-116022A
CAS No.:	333338-18-4
Molecular Formula:	C ₆ H ₁₆ NNa ₂ O ₁₂ P
Molecular Weight:	371.14
Target:	Biochemical Assay Reagents
Pathway:	Others
Storage:	-20°C, sealed storage, away from moisture * In solvent : -80°C, 6 months; -20°C, 1 month (sealed storage, away from moisture)

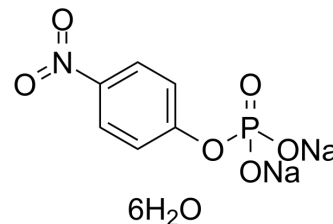

SOLVENT & SOLUBILITY

In Vitro

H₂O : 100 mg/mL (269.44 mM; Need ultrasonic)
DMSO : 10 mg/mL (26.94 mM; ultrasonic and warming and heat to 80°C)

Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		Concentration	1 mg	5 mg	10 mg
	1 mM		2.6944 mL	13.4720 mL	26.9440 mL
	5 mM		0.5389 mL	2.6944 mL	5.3888 mL
	10 mM		0.2694 mL	1.3472 mL	2.6944 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

4-Nitrophenyl phosphate (p-nitrophenyl phosphate) disodium hexahydrate is widely used as a small molecule phosphotyrosine-like substrate in activity assays for protein tyrosine phosphatases. 4-Nitrophenyl phosphate disodium hexahydrate is a colorless substrate that upon hydrolysis is converted to a yellow 4-nitrophenolate ion that can be monitored by absorbance at 405 nm^[1].

In Vitro

4-Nitrophenyl phosphate (PNPP) disodium hexahydrate is a commonly used substrate for alkaline phosphatases (ALPs). 4-Nitrophenyl phosphate disodium hexahydrate is hydrolyzed by ALP to PNP (p-nitrophenol), which quenches the fluorescence of novel gold nanoclusters (AuNCs) by the inner filter effect (IFE)^[2]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Int J Mol Sci. 2022 Feb 26;23(5):2604.

-
- J Chromatogr A. 2023 Nov 17, 464511.

See more customer validations on www.MedChemExpress.com

REFERENCES

- [1]. Lountos GT, et al. Structural analysis of human dual-specificity phosphatase 22 complexed with a phosphotyrosine-like substrate. Acta Crystallogr F Struct Biol Commun. 2015;71(Pt 2):199-205.
- [2]. Qi S, et al. Development of a facile and sensitive method for detecting alkaline phosphatase activity in serum with fluorescent gold nanoclusters based on the inner filter effect. Analyst. 2020;145(11):3871-3877.
-

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA