



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Cyclooctyne-O-PFP ester

Cat. No.:	HY-126518
CAS No.:	886209-60-5
Molecular Formula:	C ₁₆ H ₁₃ F ₅ O ₃
Molecular Weight:	348.26
Target:	ADC Linker
Pathway:	Antibody-drug Conjugate/ADC Related
Storage:	Pure form -20°C 3 years 4°C 2 years In solvent -80°C 6 months -20°C 1 month

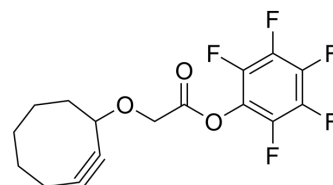

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (287.14 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM		2.8714 mL	14.3571 mL	28.7142 mL
		5 mM		0.5743 mL	2.8714 mL	5.7428 mL
		10 mM		0.2871 mL	1.4357 mL	2.8714 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (7.18 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: 2.5 mg/mL (7.18 mM); Suspended solution; Need ultrasonic					
	3. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (7.18 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	Cyclooctyne-O-PFP ester is a cleavable ADC linker used in the synthesis of antibody-drug conjugates (ADCs) ^[1] . Cyclooctyne-O-PFP ester is a click chemistry reagent, it contains an Alkyne group and can undergo copper-catalyzed azide-alkyne cycloaddition (CuAAC) with molecules containing Azide groups.
IC ₅₀ & Target	Cleavable Linker
In Vitro	ADCs are comprised of an antibody to which is attached an ADC cytotoxin through an ADC linker.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Sun, Ying, et al. Anti-PSMA Antibodies Conjugated to Nuclear Receptor Ligand Polypeptides. 20150152187A1.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA