



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

DC0-NH2

Cat. No.:	HY-129379
CAS No.:	615538-51-7
Molecular Formula:	C ₃₁ H ₂₄ ClN ₅ O ₃
Molecular Weight:	550.01
Target:	ADC Cytotoxin
Pathway:	Antibody-drug Conjugate/ADC Related
Storage:	-20°C, stored under nitrogen
	* In solvent : -80°C, 6 months; -20°C, 1 month (stored under nitrogen)

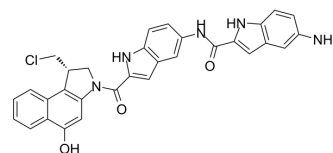

SOLVENT & SOLUBILITY

In Vitro

DMSO : 50 mg/mL (90.91 mM; Need ultrasonic)

	Solvent Concentration	Mass	1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		1.8181 mL	9.0907 mL	18.1815 mL
	5 mM		0.3636 mL	1.8181 mL	3.6363 mL
	10 mM		0.1818 mL	0.9091 mL	1.8181 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

DC0-NH2 is an effector moiety for ADC and a simplified analog of DC1 with better stability. DC0-NH2 is about 1000-fold more cytotoxic than commonly used anticancer agents (ex. Doxorubicin). DC0-NH2 can bind to the minor groove of DNA, followed by alkylation of adenine residues by its propabenzindole (CBI) component^[1].

IC₅₀ & Target

Duocarmycins

In Vitro

DC0-NH2 (0-3 nM; 72 hours) is highly potent against Ramos, Namalwa, and HL60/s cells with IC₅₀ values in the 1 pM to 10 pM range, and has 100 pM range when tested on COLO 205 cells^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Cell Viability Assay^[1]

Cell Line:	Ramos, Namalwa, HL60/s, and COLO 205 Cancer Cells
Concentration:	0-3 nM
Incubation Time:	72 hours

Result:	Inhibited Namalwa and HL60/s cells with IC ₅₀ s of 7 and 30 pM, respectively.
---------	--

REFERENCES

[1]. Zhao RY, et al. Synthesis and biological evaluation of antibody conjugates of phosphate prodrugs of cytotoxic DNA alkylators for the targeted treatment of cancer. J Med Chem. 2012 Jan 26;55(2):766-82.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA