



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Fmoc-NH-PEG8-CH2COOH

Cat. No.:	HY-133063
CAS No.:	868594-52-9
Molecular Formula:	C ₃₃ H ₄₇ NO ₁₂
Molecular Weight:	649.73
Target:	ADC Linker; PROTAC Linkers
Pathway:	Antibody-drug Conjugate/ADC Related; PROTAC
Storage:	-20°C, protect from light * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light)

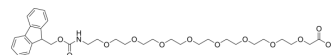

SOLVENT & SOLUBILITY

In Vitro

DMSO : 100 mg/mL (153.91 mM; Need ultrasonic)

Solvent	Mass	Concentration		
		1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM	1.5391 mL	7.6955 mL	15.3910 mL
	5 mM	0.3078 mL	1.5391 mL	3.0782 mL
	10 mM	0.1539 mL	0.7696 mL	1.5391 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

Fmoc-NH-PEG8-CH2COOH is a cleavable ADC linker used in the synthesis of antibody-drug conjugates (ADCs). Fmoc-NH-PEG8-CH2COOH is also a PEG-based PROTAC linker that can be used in the synthesis of PROTACs^[1].

IC₅₀ & Target

Cleavable Linker

PEGs

Alkyl/ether

In Vitro

ADCs are comprised of an antibody to which is attached an ADC cytotoxin through an ADC linker. PROTACs contain two different ligands connected by a linker; one is a ligand for an E3 ubiquitin ligase and the other is for the target protein. PROTACs exploit the intracellular ubiquitin-proteasome system to selectively degrade target proteins. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Michael A, et al. Synthesis of Bifunctional Integrin-Binding Peptides Containing PEG Spacers of Defined Length for Non-Viral Gene Delivery. Volume2008, Issue17.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA