



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com


www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



IR-820

Cat. No.:	HY-136886
CAS No.:	172616-80-7
Molecular Formula:	$C_{46}H_{50}ClN_2NaO_6S_2$
Molecular Weight:	849.47
Target:	Fluorescent Dye
Pathway:	Others
Storage:	4°C, sealed storage, away from moisture * In solvent : -80°C, 6 months; -20°C, 1 month (sealed storage, away from moisture)

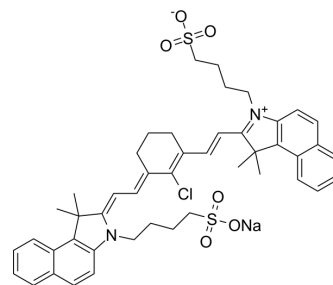

SOLVENT & SOLUBILITY

In Vitro

DMSO : 5 mg/mL (5.89 mM; ultrasonic and warming and heat to 60°C)

	Solvent Concentration	Mass	1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		1.1772 mL	5.8860 mL	11.7720 mL
	5 mM		0.2354 mL	1.1772 mL	2.3544 mL
	10 mM		---	---	---

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

IR-820 (New Indocyanine Green) is an infrared blood pool contrast agent. IR-820 also is normally used as a laser and near-infrared dye to detect and quantify diseased tissue in live animals^[1].

In Vitro

Guidelines (Following is our recommended protocol. This protocol only provides a guideline, and should be modified according to your specific needs).

In vivo^[1]:

1. The powdered IR-820 was freshly mixed with phosphate buffer saline (PBS) to a final concentration of 0.2 mM and 100 µl of this solution was used for injecting into experimental animals.
2. Used 2 groups of four Hairless SKHi/sKHL mice. The first group was injected with IR-820 by intravenous tail vein injection while the second was injected intraperitoneally.
3. IR-820 has maximal excitation and emission wavelengths of 710 nm and 820 nm, respectively. Both groups were serially imaged for 8 days.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

-
- Int J Nanomedicine. 2023 Oct 31.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Suresh I Prajapati, et al. Near-infrared imaging of injured tissue in living subjects using IR-820. Mol Imaging. 2009 Jan-Feb;8(1):45-54.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA