



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)

Levamlodipine-d₄

Cat. No.:	HY-14744S
CAS No.:	1346617-19-3
Molecular Formula:	C ₂₀ H ₂₁ D ₄ ClN ₂ O ₅
Molecular Weight:	412.9
Target:	Calcium Channel; Isotope-Labeled Compounds
Pathway:	Membrane Transporter/Ion Channel; Neuronal Signaling; Others
Storage:	-20°C, protect from light, stored under nitrogen * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light, stored under nitrogen)

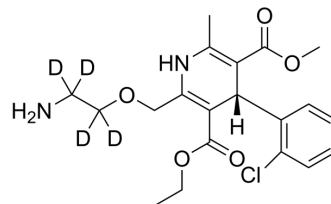

SOLVENT & SOLUBILITY

In Vitro

DMSO : 100 mg/mL (242.19 mM; Need ultrasonic and warming)

Concentration	Mass		
	1 mg	5 mg	10 mg
1 mM	2.4219 mL	12.1095 mL	24.2189 mL
5 mM	0.4844 mL	2.4219 mL	4.8438 mL
10 mM	0.2422 mL	1.2109 mL	2.4219 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

Levamlodipine-d₄ is the deuterium labeled Levamlodipine. Levamlodipine ((S)-Amlodipine) is a powerful dihydropyridine calcium channel blocker, possessing vasodilation properties and used in the treatment of hypertension and angina[1].

In Vitro

Stable heavy isotopes of hydrogen, carbon, and other elements have been incorporated into drug molecules, largely as tracers for quantitation during the drug development process. Deuteration has gained attention because of its potential to affect the pharmacokinetic and metabolic profiles of drugs^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Russak EM, et al. Impact of Deuterium Substitution on the Pharmacokinetics of Pharmaceuticals. *Ann Pharmacother.* 2019;53(2):211-216.

[2]. Jain R, et al. Polyaniline-graphene oxide nanocomposite sensor for quantification of calcium channel blocker levamlodipine. *Mater Sci Eng C Mater Biol Appl.* 2016 Aug 1;65:205-14.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA