



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Hoechst S 769121

Cat. No.:	HY-15619
CAS No.:	74681-68-8
Molecular Formula:	C ₂₅ H ₂₈ Cl ₃ N ₇ O ₂ S
Molecular Weight:	596.96
Target:	Fluorescent Dye
Pathway:	Others
Storage:	4°C, sealed storage, away from moisture and light * In solvent : -80°C, 6 months; -20°C, 1 month (sealed storage, away from moisture and light)

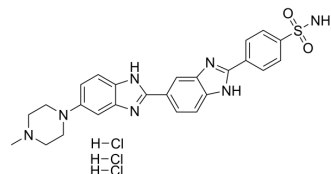

SOLVENT & SOLUBILITY

In Vitro

H₂O : ≥ 6 mg/mL (10.05 mM)
* "≥" means soluble, but saturation unknown.

Concentration	Solvent	Mass		
		1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM	1.6752 mL	8.3758 mL	16.7515 mL
	5 mM	0.3350 mL	1.6752 mL	3.3503 mL
	10 mM	0.1675 mL	0.8376 mL	1.6752 mL

Please refer to the solubility information to select the appropriate solvent.

In Vivo

1. Add each solvent one by one: PBS
Solubility: 5 mg/mL (8.38 mM); Clear solution; Need ultrasonic and warming and heat to 60°C

BIOLOGICAL ACTIVITY

Description

Hoechst S 769121 is a marker dye in Hoechst series. Hoechst is a live nuclear marker dye. Hoechst binds to the grooves in the DNA double strand, which tends to be A/T-rich DNA strand. Although it binds to all nucleic acids, the A/T-rich double strand DNA significantly enhances fluorescence intensity. Therefore, Hoechst dye can be used for living cell labeling. The fluorescence intensity of Hoechst dye increases with the increase of pH of solution^[1].

In Vitro

General Protocol
Preparation of Hoechst working solution
1.1 Preparation of the stock solution
Dissolve 10 mg of in 5 mL DMSO
Note: It is recommended to store the stock solution at 4°C or -20°C away from light and avoid repetitive freeze-thaw cycles.
1.2 Preparation of Hoechst working solution
Dilute the stock solution in serum-free cell culture medium or PBS to obtain final concentration 10 µg/mL Hoechst working

solution.

Note: Please adjust the concentration of Hoechst working solution according to the actual situation.

1. Cell staining

2.1 Suspension cells 6-well plate

a. Centrifuge at 1000 g at 4°C for 3-5 minutes and then discard the supernatant. Wash twice with PBS, 5 minutes each time.

The cell density is 1×10^6 /mL.

b. Add 1 mL of working solution, and then incubate at room temperature for 3-10 minutes.

c. Centrifuge at 400 g at 4°C for 3-4 minutes and then discard the supernatant.

d. Wash twice with PBS, 5 minutes each time.

e. Resuspend cells with serum-free cell culture medium or PBS. Observation by fluorescence microscopy or flow cytometry.

2.2 Adherent cells

a. Culture adherent cells on sterile coverslips.

b. Remove the coverslip from the medium and aspirate excess medium.

c. Add 100 µL of working solution, gently shake it to completely cover the cells, and then incubate at room temperature for 3-10 minutes.

d. Wash twice with medium, 5 minutes each time. Observation by fluorescence microscopy or flow cytometry.

Precautions

1. Please adjust the concentration of Hoechst working solution according to the actual situation.

2. This product is for R&D use only, not for drug, household, or other uses.

3. For your safety and health, please wear a lab coat and disposable gloves to operate.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Latt SA, Stetten G, Juergens LA, Recent developments in the detection of deoxyribonucleic acid synthesis by 33258 Hoechst fluorescence. The journal of histochemistry and cytochemistry : official journal of the Histochemistry Society 23 (7): 493-505.

[2]. a b c "Hoechst Stains". Invitrogen (Molecular Probes).

[3]. Portugal J, Waring MJ. Assignment of DNA binding sites for 4',6-diamidino-2-phenylindole and bisbenzimidazole (Hoechst 33258). A comparative footprinting study. Biochimica et Biophysica Acta 949

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA