



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Antibacterial agent 206

Cat. No.:	HY-162430
Molecular Formula:	C <sub>30</sub> H <sub>28</sub> FN <sub>5</sub> O <sub>5</sub>
Molecular Weight:	557.57
Target:	Bacterial
Pathway:	Anti-infection
Storage:	Please store the product under the recommended conditions in the Certificate of Analysis.

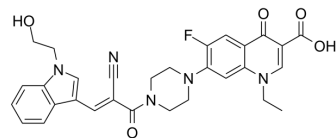

### BIOLOGICAL ACTIVITY

<b>Description</b>	Antibacterial agent 206 (Compound 10e) is a indolylacryloyl-derived oxacins, which exhibits broad antibacterial spectrum with MIC of 0.25-1 µg/mL. Antibacterial agent 206 reduces the exopolysaccharide, eliminates the biofilm, and thus attenuates the drug resistance. Antibacterial agent 206 exhibits antibacterial activity through destory of membrane integrity, accumulation of reactive oxygen species ROS, and inhibition of DNA replication <sup>[1]</sup> .								
<b>In Vitro</b>	<p>Antibacterial agent 206 exhibits antimicrobial activity against Enterococcus, Staphylococcus, Escherichia, Acinetobacter and Pseudomonas strains with MIC between 0.25-1 µg/mL<sup>[1]</sup>.</p> <p>Antibacterial agent 206 (0-256 µg/mL) exhibits low hemolysis, low cytotoxicity in blood, human hepatocytes (LO2) and HUVECs, and good compatibility with red blood cells<sup>[1]</sup>.</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p> <p>Cell Viability Assay<sup>[1]</sup></p> <table border="1" style="width: 100%; border-collapse: collapse;"> <tr> <td style="width: 30%;">Cell Line:</td> <td>LO2, HUVEC</td> </tr> <tr> <td>Concentration:</td> <td>2-64 µg/mL</td> </tr> <tr> <td>Incubation Time:</td> <td></td> </tr> <tr> <td>Result:</td> <td>Inhibited cell viability slightly in dose-dependent manner.</td> </tr> </table>	Cell Line:	LO2, HUVEC	Concentration:	2-64 µg/mL	Incubation Time:		Result:	Inhibited cell viability slightly in dose-dependent manner.
Cell Line:	LO2, HUVEC								
Concentration:	2-64 µg/mL								
Incubation Time:									
Result:	Inhibited cell viability slightly in dose-dependent manner.								

### REFERENCES

[1]. Hu YG, et al., Discovery of indolylacryloyl-derived oxacins as novel potential broad-spectrum antibacterial candidates. Eur J Med Chem. 2024 Apr 15;270:116392.

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA