



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Biotin NHS

|                    |  |
|--------------------|--|
| Cat. No.:          | HY-D0802   |
| CAS No.:           | 35013-72-0   |
| Molecular Formula: | C <sub>14</sub> H <sub>19</sub> N <sub>3</sub> O <sub>5</sub> S  |
| Molecular Weight:  | 341.38   |
| Target:            | Fluorescent Dye  |
| Pathway:           | Others   |
| Storage:           | -20°C, sealed storage, away from moisture and light<br>* The compound is unstable in solutions, freshly prepared is recommended. |

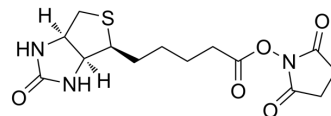

### SOLVENT & SOLUBILITY

|   |   |                          |           |            |            |
|---|---|--------------------------|-----------|------------|------------|
| In Vitro  | DMSO : 62.5 mg/mL (183.08 mM; Need ultrasonic)  |                          |           |            |            |
|   |   | Solvent<br>Concentration | Mass      |            |            |
|   | Preparing<br>Stock Solutions  | 1 mM                     | 1 mg      | 5 mg       | 10 mg      |
|   |   | 5 mM                     | 2.9293 mL | 14.6464 mL | 29.2929 mL |
|   |   | 10 mM                    | 0.5859 mL | 2.9293 mL  | 5.8586 mL  |
|   | 10 mM   | 0.2929 mL                | 1.4646 mL | 2.9293 mL  |            |
| Please refer to the solubility information to select the appropriate solvent. |   |                          |           |            |            |
| In Vivo   | 1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline<br>Solubility: ≥ 2.08 mg/mL (6.09 mM); Clear solution |                          |           |            |            |
|   | 2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)<br>Solubility: ≥ 2.08 mg/mL (6.09 mM); Clear solution            |                          |           |            |            |
|   | 3. Add each solvent one by one: 10% DMSO >> 90% corn oil<br>Solubility: ≥ 2.08 mg/mL (6.09 mM); Clear solution                            |                          |           |            |            |

### BIOLOGICAL ACTIVITY

|             |   |
|-------------|---|
| Description | Biotin NHS is an amino reactive biotin reagent used in the preparation of biotinylated surfaces or polypeptides.  |
| In Vitro    | <p>Guidelines (Following is our recommended protocol. This protocol only provides a guideline, and should be modified according to your specific needs).</p> <p>?Calculations:</p> <p>?For best results, use ≥ 12-fold molar excess of biotin for a 10 mg/mL protein solution or ≥ 20-fold molar excess of biotin for a 2 mg/mL protein solution.</p> <p>?1. Calculate millimoles of biotin reagent to add to the reaction for a 20-fold molar excess:</p> <p>?mL protein × (mg protein/mL protein) × (mmol protein/mg protein) × (20 mmol protein/mmol protein) = mmol Biotin.</p> |

?\* 20 = Recommended molar fold excess of biotin for 2 mg/ml protein sample  
?2. Calculate microliters of 10 mM biotin reagent solution (prepared in Step B.3) to add to the reaction:  
$$?mmol \text{ Biotin} \times (1000000 \mu\text{L/L}) \times (\text{L} / 10 \text{ mmol}) = \mu\text{L Biotin}$$
  
?Biotin Labeling Reaction:  
?1. Dissolve 1-10 mg protein in 0.5-2.0 mL PBS according to the calculation made in Calculations.  
?2. Before use, prepare a 10 mM solution of the biotin reagent in an organic solvent such as DMSO.  
?3. Add the appropriate volume ( $\mu\text{L}$ ) of 10 mM biotin reagent solution to the protein solution.  
?4. Incubate reaction on ice for two hours or at room temperature for 30 minutes.  
?5. To tests the labeled protein by ELISA or Western blot.  
MCE has not independently confirmed the accuracy of these methods. They are for reference only.

## CUSTOMER VALIDATION

- ACS Nano. 2023 Feb 14;17(3):2537-2553.
- ACS Nano. 2019 Feb 26;13(2):2501-2510.
- Biomacromolecules. 2019 Jan 14;20(1):469-477.
- J Virol. 2021 May 26;JVI0018721.
- Vet Microbiol. 2023 Apr 11;281:109743.

See more customer validations on [www.MedChemExpress.com](http://www.MedChemExpress.com)

## REFERENCES

- [1]. Wojchowski DM et al. Hybridoma production by simplified avidin-mediated electrofusion. J Immunol Methods. 1986 Jun 24;90(2):173-7.
- [2]. Ngo TT et al. Biotinyl-glucose-6-phosphate dehydrogenase preparation, kinetics, and modulation by avidin. Appl Biochem Biotechnol. 1982 Nov;7(6):443-54.

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA