



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

H2DCFDA

Cat. No.:	HY-D0940
CAS No.:	4091-99-0
Molecular Formula:	C ₂₄ H ₁₆ Cl ₂ O ₇
Molecular Weight:	487.29
Target:	Reactive Oxygen Species
Pathway:	Immunology/Inflammation; Metabolic Enzyme/Protease; NF-κB
Storage:	-20°C, protect from light * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light)

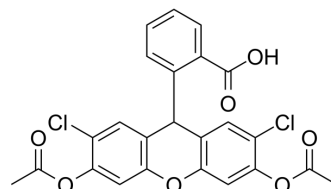

SOLVENT & SOLUBILITY

In Vitro

DMSO : 125 mg/mL (256.52 mM; Need ultrasonic)
Ethanol : 20 mg/mL (41.04 mM; Need ultrasonic)

Preparing Stock Solutions	Solvent Concentration	Mass		
		1 mg	5 mg	10 mg
	1 mM	2.0522 mL	10.2608 mL	20.5217 mL
	5 mM	0.4104 mL	2.0522 mL	4.1043 mL
	10 mM	0.2052 mL	1.0261 mL	2.0522 mL

Please refer to the solubility information to select the appropriate solvent.

In Vivo

- Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline
Solubility: 2.08 mg/mL (4.27 mM); Suspended solution; Need ultrasonic
- Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)
Solubility: ≥ 2.08 mg/mL (4.27 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% corn oil
Solubility: ≥ 2.08 mg/mL (4.27 mM); Clear solution

BIOLOGICAL ACTIVITY

Description

H2DCFDA (DCFH-DA) is a cell-permeable probe used to detect intracellular reactive oxygen species (ROS) (Ex/Em=488/525 nm)^[1].

In Vitro

Guidelines (Following is our recommended protocol. This protocol only provides a guideline, and should be modified according to your specific needs).

- H2DCFDA is dissolved in DMSO to obtain a 10 mM stock solutions and further diluted before use.
- Cells are incubated with 5 μM H2DCFDA solution in PBS in the dark for 30 min at 37°C, then harvested with 0.05% trypsin-EDTA solution, suspended in a fresh medium, and immediately analyzed with flow cytometer.

3. Along with the H2DCFDA probe, if indicated, use ROS-insensitive modification of the fluorescein dye DCFDA as a positive control. The staining procedure is the same as for the H2DCFDA^[1].
MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Kinase Assay

ROS Measurements^[1]

For the detection of intracellular ROS level, ROS-sensitive probe H2DCFDA is used. Adherent cells (ESCs, difESCs, eMSCs, HeLa, U118) are incubated with 5 μ M staining solution in PBS in the dark for 30 min at 37°C, then harvested with 0.05% trypsin-EDTA solution, suspended in a fresh medium, and immediately analyzed with flow cytometer. Lymphocytes, both control and PHA-activated, are resuspended in PBS, incubated with 5 μ M of H2DCFDA in the dark for 30 min at 37°C, and immediately analyzed. Along with the H2DCFDA probe, if indicated, ROS-insensitive modification of the fluorescent dye DCFDA is used. The staining procedure is the same as for the H2DCFDA^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Cell. 2021 Jun 24;184(13):3528-3541.e12.
- Signal Transduct Target Ther. 2022 Aug 15;7(1):288.
- Adv Mater. 2023 Sep 5;e2306469.
- Drug Resist Updat. 2023 Jul;69:100974.
- Cell Mol Immunol. 2023 Aug 9.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Lyublinskaya OG, et al. Redox environment in stem and differentiated cells: A quantitative approach. Redox Biol. 2017 Aug;12:758-769.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA