



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

5-DTAF

Cat. No.:	HY-D1082
CAS No.:	51306-35-5
Molecular Formula:	C ₂₃ H ₁₂ Cl ₂ N ₄ O ₅
Molecular Weight:	495.27
Target:	Fluorescent Dye
Pathway:	Others
Storage:	-20°C, protect from light * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light)

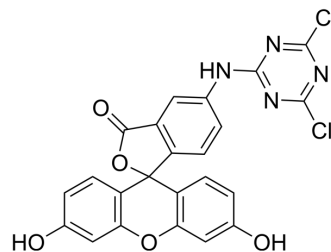

SOLVENT & SOLUBILITY

In Vitro

DMSO : 50 mg/mL (100.96 mM; Need ultrasonic)

Concentration	Solvent	Mass		
		1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM	2.0191 mL	10.0955 mL	20.1910 mL
	5 mM	0.4038 mL	2.0191 mL	4.0382 mL
	10 mM	0.2019 mL	1.0096 mL	2.0191 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description

5-DTAF is a isomeric fluorescein derivative that directly react with polysaccharides and other alcohols in aqueous solution at pH 9. The maximum of the 5-DTAF emission profile is at a wavelength of 518 nm after excitation at 488 nm^[1].

REFERENCES

[1]. Fergg F, et al. Investigations of the microscopic structure of poly (vinyl alcohol) hydrogels by confocal laser scanning microscopy. Colloid and Polymer Science, 2001, 279(1): 61-67.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA