



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Hosenkoside K

Cat. No.:	HY-N2243		
CAS No.:	160896-49-1		
Molecular Formula:	C ₅₄ H ₉₂ O ₂₅		
Molecular Weight:	1141.29		
Target:	Others		
Pathway:	Others		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	6 months
		-20°C	1 month

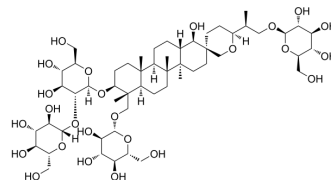

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (87.62 mM; Need ultrasonic)			
		Solvent Concentration	Mass	
			1 mg	5 mg
			10 mg	
Preparing Stock Solutions	1 mM	0.8762 mL	4.3810 mL	8.7620 mL
	5 mM	0.1752 mL	0.8762 mL	1.7524 mL
	10 mM	0.0876 mL	0.4381 mL	0.8762 mL
Please refer to the solubility information to select the appropriate solvent.				
In Vivo	<ol style="list-style-type: none"> Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (2.19 mM); Clear solution Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.5 mg/mL (2.19 mM); Clear solution Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (2.19 mM); Clear solution 			

BIOLOGICAL ACTIVITY

Description	Hosenkoside K is a baccharane glycoside isolated from the seeds of <i>impatiens balsamina</i> .
In Vivo	The mean peak concentrations of hosenkoside K in rats after oral administration of total saponins of <i>Semen Impatiensis</i> is 511.11±234.07 ng/mL, occurring at times of 0.46 h. The terminal half-lives of hosenkoside A is 4.96±1.75 h ^[1] . MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Animal Administration

Male Wistar rats (7 weeks old, 220±20 g, n=6) are fasted for 12 h. Total saponins of Semen Impatientis at a dose of 300 mg (in 10 mL)/kg are orally administered to rats using a gastric gavage tube. The doses of hosenkoside A is 12.4 mg/kg for an oral dose of 300 mg/kg total saponins. Blood samples (0.20 mL) are collected from suborbital vein into heparinized tubes before administration (0 min), and 0.08, 0.25, 0.5, 0.75, 1, 1.5, 2, 4, 6, 8, 12 and 24 h after dosing.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Yu X, et al. An LC-MS/MS method for simultaneous determination of hosenkoside A and hosenkoside K from Semen Impatientis in rat plasma and its application to a pharmacokinetic study. *Biomed Chromatogr.* 2017 Mar;31(3).

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA